




Diagnosis

Explorer	2
Excavator	5
Dressing Tweezers	6
Mouth Mirror	8
Disposable Instrument	10
Temporary Filling Instrument	10
Evacuation Tip	11
Stainless Dish / Tray	12



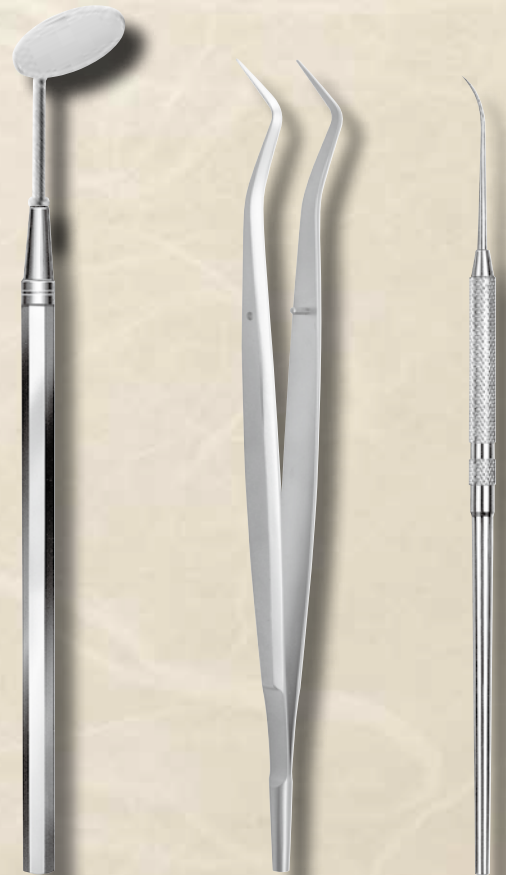

Explorer

-Single End-



#3
No. 01003



#5
No. 01005




#9
No. 01009



#23
No. 01023



#25
No. 01025

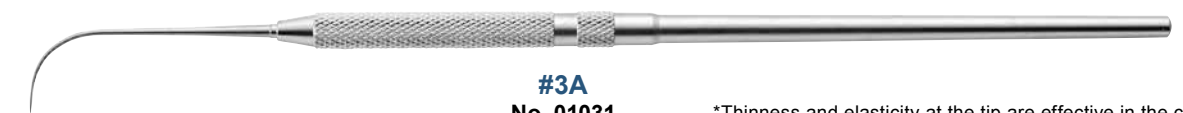

#3A
No. 01031

*Thinness and elasticity at the tip are effective in the check of dental tartar, cavity and surplus cement.



Multi
No. 01030


*A tip part is made thick for intensity.



#23BS
No. 01024



Octagonal Handle #9
No. 01109

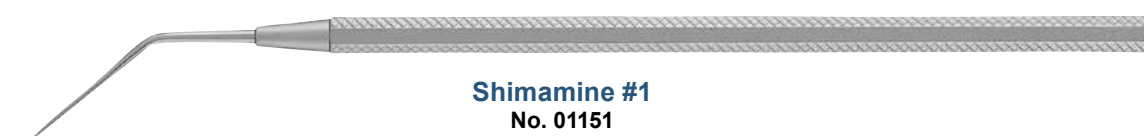


Shimamine #1
No. 01151

-Ball Point-




#1
No. 01221

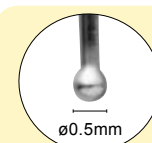


*The round tip realizes checking of caries without damaging pit and fissure tissues.
*Suitable for filling of sealant.


-Double End-

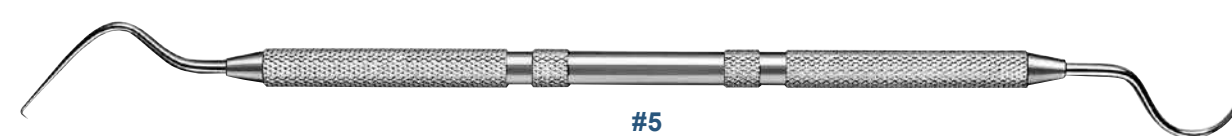

#5
No. 01205

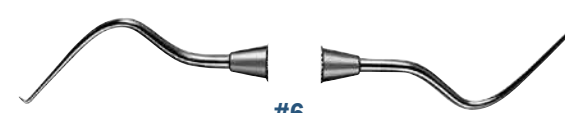

#6
No. 01206

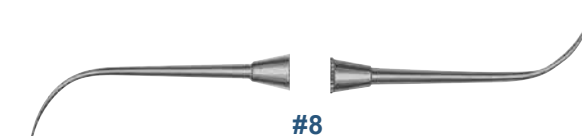

#8
No. 01208



#3-9
No. 01209



#11-12
No. 09224



*Explorer with the shank of "Gracey Curette #11-12".
*Thinness and elasticity at the tip are effective in the check of dental tartar, cavity and surplus cement.

- for Root Canal Therapy -




#16
No. 01316





Dressing Tweezers

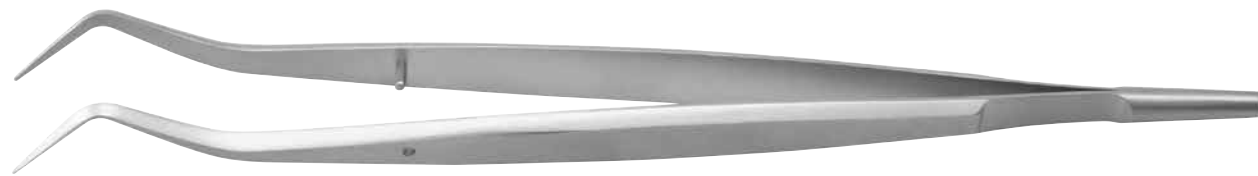


Koshi
No. 11020

Size / 165mm



NC #18 with Diamond
No. 11307

Size / 160mm

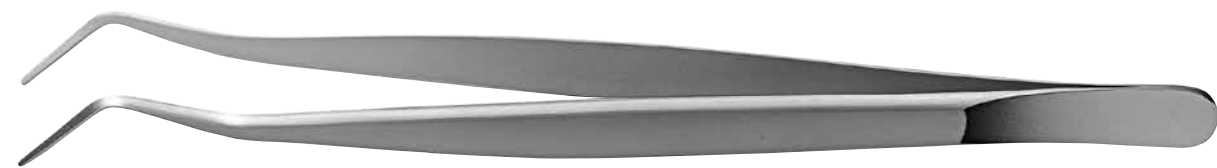

Premium
No. 11019


Size / 160mm

*Delicate tweezers.
*Careful finish prevents tip misalignment.




Mini #18
No. 11308

Size / 145mm

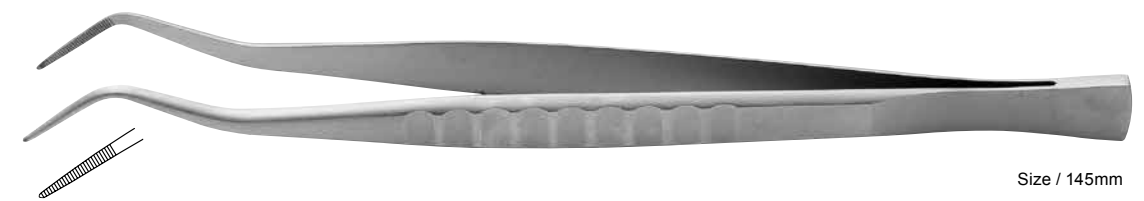

Mini #18 with Serration
No. 11309

Size / 145mm

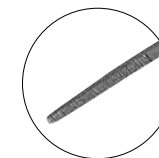

*Serrated tip for stable hold.



Size / 126mm

Multi Use Tweezers Fine
No. 27409

**"Multi" use for handling bracket, tube, lingual button and so on.

*Firm and stable hold without dropping off even used for tube thanks to its "exquisite" tip design.

*Designed slim tip not to get stuck with the space of hook part of tube.


[Handling bracket]



[Handling tube]



*Designed angle for clear view.

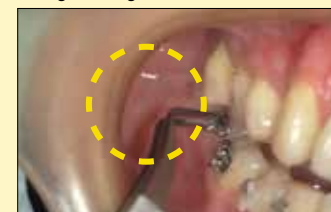

*Bonding tube on second molar.
Designed tip angle not to contact retractor.

*Compact design even for lingual treatment..




*Able to use for holding and tying elastic tube.



-Clear 2-

*With thin glass, there is little blur of an image and reflected vividly.

*Suitable only for YDM mirror handle.

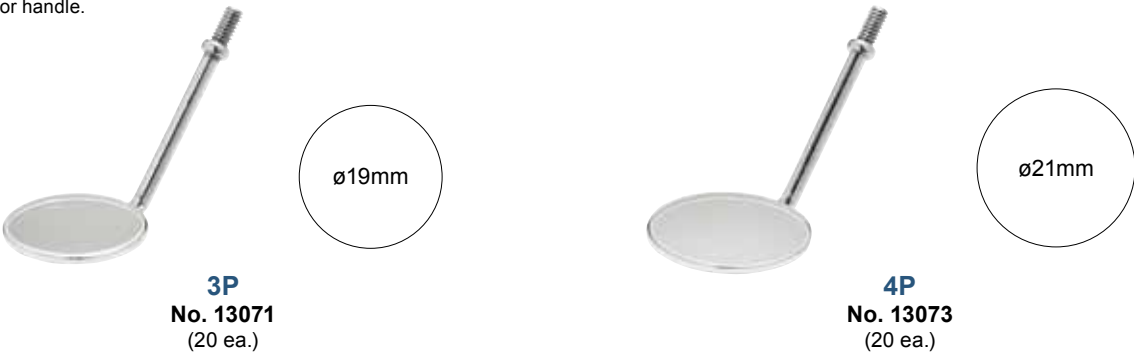

*Thickness / 1.9mm


*Reflection rate / About 93%

-Front Surface-

***Titanium Coating**

*Direct reflection prevents an image blurring.

*Suitable only for YDM mirror handle.

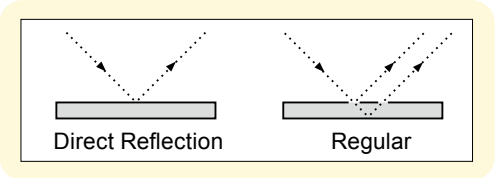


*Thickness / 1.9mm

*Reflection rate / About 98%

Metal Mirror (Simple Stem)

*The mirror made from stainless steel does not have a fear of being divided, and is suitable for the medical treatment of child dentistry, disabled person dentistry.

*Suitable only for YDM mirror handle.

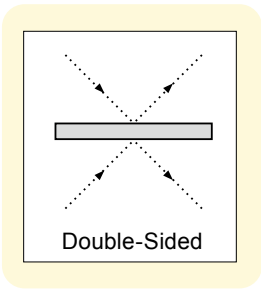

No. 13131
(6 ea.)




*Thickness / 1.1mm

*Reflection rate / About 65%

Mirror Handle

*Suitable only for YDM mouth mirror.

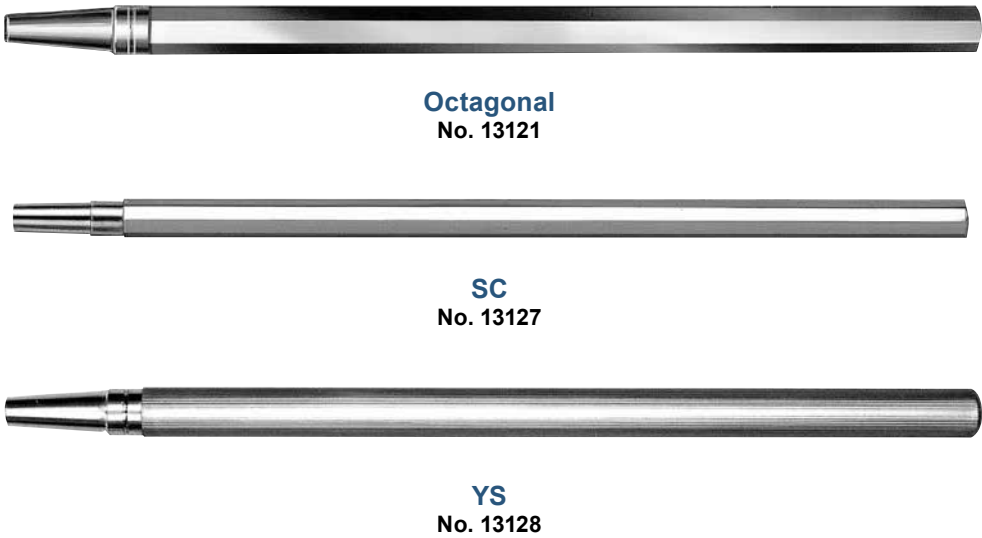

Mirror Handle with Millimeter Scale

*Suitable only for YDM mouth mirror.



*Material / POM

[How to measure endodontic file or reamer.]



①Fit stopper of reamer or file and pull in the direction of arrow.



②Measure the length of reamer or file with scale on the handle.



③Stopper of reamer or file is fixed.


*Measurement: Up to 30mm



Mouth Mirror
No. 28102
(50 ea.)



Explorer
No. 28101
(50 ea.)

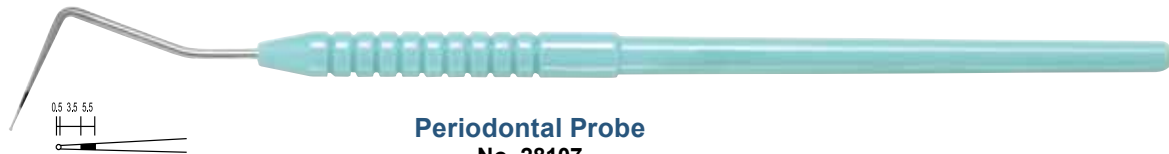

Periodontal Probe
No. 28107
(50 ea.)

*Disposable type (Not sterilized)
*Material / PS : Heatproof temperature 70°C (Handle)
*Need disinfection by solution or EOG sterilization before use.



Temporary Filling Carrier

No. 13200



*A tip is heated and used.
*To prevent the burn of membrane, a black cover is made to slide to a tip.



Temporary Filling Remover




#10
No. 03010

Evacuation Tip

-Universal-

*There is no trap for suction prevention attached at the evacuation tip, so it is easy to clean. (Rubber sleeve has a trap)

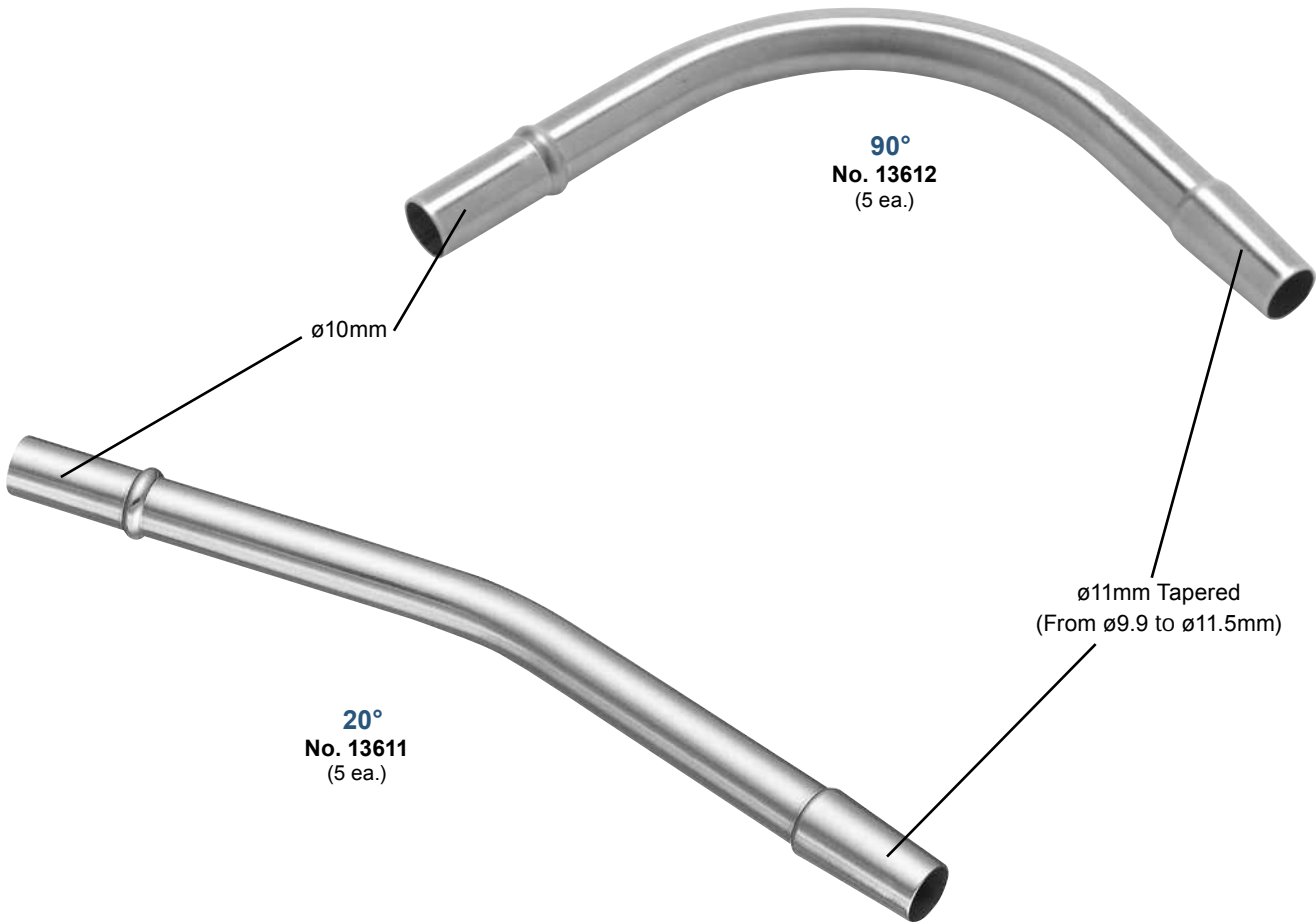

Sold Separately

Rubber Sleeve
No. 13613
(5 ea.)



*Size / 40mm
*Material / Si



*"Rubber Sleeve" has a trap inside.

Adapter
No. 13614



*Material / PTFE

*In case the inside diameter of evacuation tube is 16mm, this adaptor is needed for the use of "Evacuation Tip Universal".

Cleaning Brush L (12mm)
No. 13653
(12 ea.)



Size / 180mm



Stainless Dish



*Excellent when using small portion of bath such as topical anesthesia, local disinfection, etching solution, and temporary filling material (powder or liquid).

No. 13670
(5 ea.)



*Hole diameter / $\phi 16\text{mm}$
*Size / W56 x D31 x H8.5mm

Stainless Tray



M
No. 13413



Size / W245 x D146 x H8mm

S
No. 13410



Size / W245 x D73 x H8mm

