


Endodontics

Rubber Dam Product	52
Endodontic Product	54
Fiber Cutter	62
GP Remover Spear	63

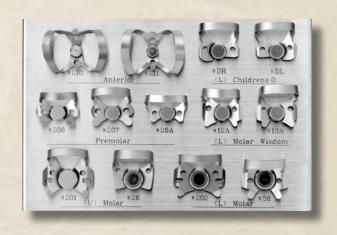

Rubber Dam Clamp

*Stainless Steel

*#12A / 13A \rightarrow Serrated jaws and offset jaw height provides stability to partially erupted teeth.

*#26 / #26A → Used in pedodontic. In that case, #26 is for deciduous molars and #26A is for deciduous anterior teeth.



#210 Anterior teeth No. 08510



#211 Anterior teeth No. 08511



Lower first deciduous molar right No. 08560



Lower first deciduous molar left No. 08561



#206 Premolars No. 08506



Premolars No. 08507



#26A Premolars No. 08527



No. 08512



Lower molars / Wisdom teeth Lower molars / Wisdom teeth No. 08513



#201 Upper molars No. 08501



#26 Upper molars No. 08526



#202 Lower molars No. 08502



#56 Lower molars No. 08556

-Clamp shape-



Fit on gingival margin.



Gum is not pressed. (Grasp part is flat)

Board No. 08258



- How to use -

Size / W132 x D87mm

Young Frame








Reamer Guard





- *Can be arranged "Reamer Guard" by reamer or file type and size.
- *Cover will not allow the reamers and files to fall out.
- *An Engine reamer can be stored with "Bur Guard".

No.	Description
08431	Yellow
08432	Orange
08433	Pink
08434	Blue
08435	Green

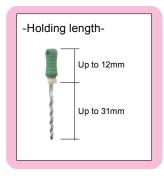

- *Storage capacity / 8 reamers and files
- *Size / W74 x D10 x H52 mm

- *Sterilizable with liquid sterilizing solutions.
 *Product does not include reamers and files.







*Has a built in gauge to measure stop positions during root canal treatments.




*Able to connect each "Reamer Guards".

Reamer Stand

-Safety-

- *It is ideal for managing reamers and files for every patients.
- *Can be autoclave or ultrasonic cleaning with reamers and files.
- *Lid does not allow reamers and files to fall out.

No. 08415

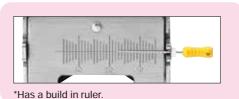

- *Storage capacity / 11 reamers and files *Holding length / From 21 to 31mm *Size / W48 x D43 x H47mm

- *Material / PP (Green part)
 *Product does not include reamers and files.



Clean Pot



Contents Main Body / Paper (100 sheets)

No. 08423

- *During the treatment, can be cleaned reamers and files by paper.
- *Set 2 sheets to stabilize reamers for cleaning.



Paper No. 08422 (300 sheets)



55

^{*}Size / ø45 x H36mm

^{*}Product does not include reamers and files.

Endodontic Aspirator Type-2



Contents Main Body / Joint (ø10mm & ø11mm Tapered) / Nozzle (5 ea.) / Cleaning Wire



*The inside of a root canal is dried by an air current with suction in a root canal.

*Because of drying by an air current, it also becomes saving of a paper point.



*Tip outside diameter / ø0.6mm









*In case the inside diameter of evacuation tube is 16mm, this adaptor is needed for the use of Evacuation Tip Universal.

Irrigating Syringe







Explorer







for Root Canal Therapy #16 No. 01316

Broach Holder





File Holder Type-2



- *Excellent for use in hard-to-reach posterior teeth.
- *Applicable for root canal treatment of both mesial and distal.







Mesial

Distal

Finger Ruler



Scale



No. 21155





Endodontic Gauge



57

No. 21154 Front side




X-Ray Gauge

Contents Sleeve (2 ea.) / Gauge (2 ea.)

No. 21157





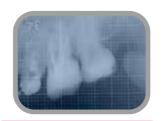

Gauge (1mm by 1mm with hold line every 5mm) *Suitable for root canal treatment, accurate examination before placing an implant, checking of calculus quantity and its area.



Gutta Percha Case

*Size / W38 x D48 x H8.5mm (Sleeve) W33 x D45mm (Gauge) *Material / Si (Sleeve)





*Can be autoclave with points.

Mini No. 24017


*Size / W91 x D42 x H13mm

*Material / Si (Green part)
*Product does not include main and accessory points.



Endodontic Tweezers




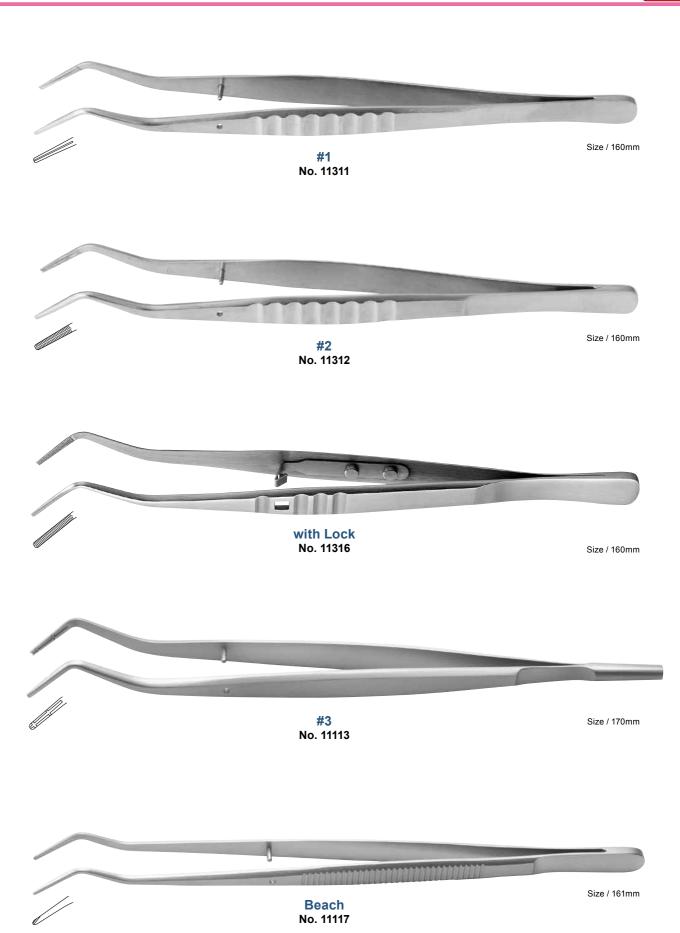



Spreader



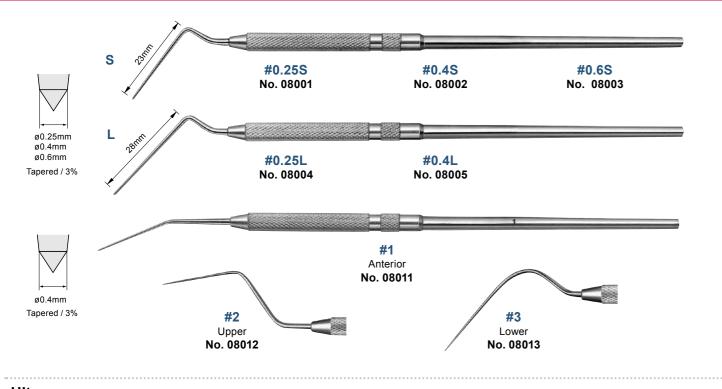

-Ultra*Superelastic Alloy


*Since an elastic metal is used, the root canal which is curving can also be pressurized in accordance with a root wall.



Spreader Tip (Replaceable Type)




-Handle-

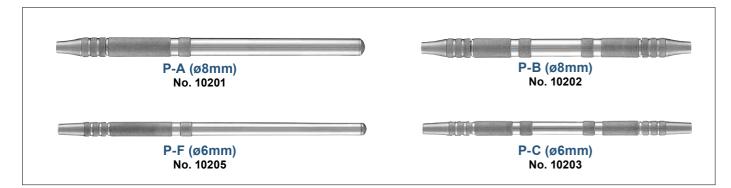





Root Canal Plugger



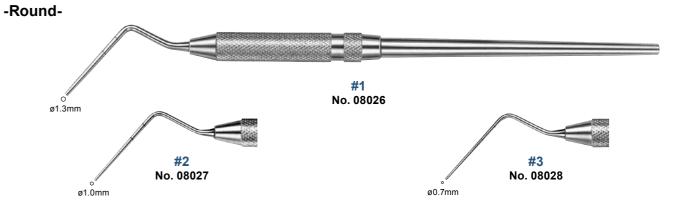

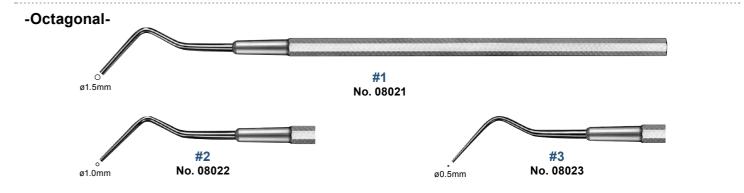

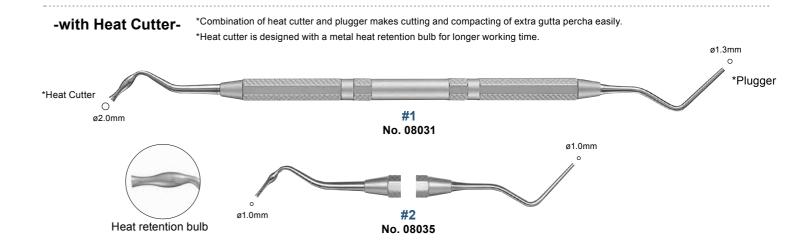

Heat Cutter



Fiber Cutter Type2

Contents Main body / Holder / Sleeve Guide (3 ea.)

- *Easy to cut the fiber post length for your need.
- *There is no fluff of a cutting plane.
- *Particulates are not generated.

No. 08071



- *Cut fiber can be entered in to the
- *Easy to cut the fiber thanks to the gauge on the holder.
- (Measurement from 5mm to 15mm)



Fiber Cutter Type2	Nipper	Steal Round Bur	Cutting Disc	Cutting Disc with Diamond

Problems from using Steal Round Bur, Cutting Disc, Diamond Disc

- 1) Fiber become fluffy.
- 2) Melt the surface because of the friction heat between fiber post and bur.
- 3) Remain the color of bur or disc on the section.

GP Remover Spear



*Remove all joint parts and replaceable tip before sterilization and disinfection.

*For more replaceable type products see P. 211.




*In Re-treatment of root canal, old Gutta Percha must be removed.


*After using ultrasonic, old Gutta Percha can be soften and then it can be removed by GP Remover Spear.



Tips' shape is tiny and have 360dergees of edges, so easy to achieve the target point. It can be used as finalizing remover after using rotary files.




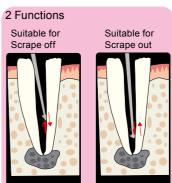

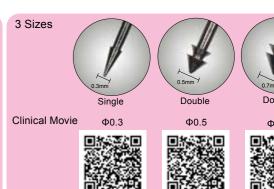






No.	standard	Tip diameter×Length	Tips shape
27131	0.38	Ф0.3×21mm	Single
27132	0.3L	Ф0.3×28mm	
27133	0.5S	Ф0.5×21mm	
27134	0.5L	Ф0.5×28mm	Double
27135	0.7S	Ф0.7×21mm	40
27136	0.7L	Ф0.7×28mm	

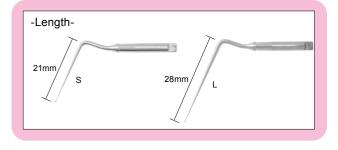

-Handle-

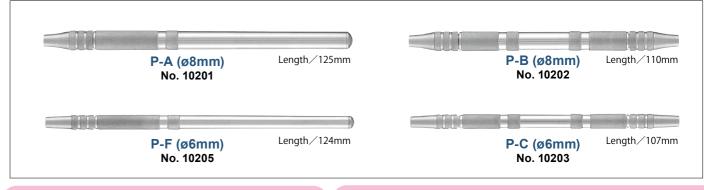

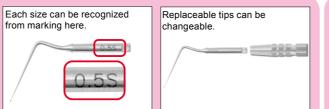

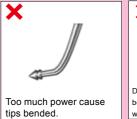





^{*}Size / 130mm

^{*}Material / TPX (Holder)

