




Surgery

Hemostatic Forceps	66
Surgical Curette	67
Elevator	68
Periosteal Elevator	76
Periodontal Knife	81
Mallet	82
Retractor	83
Gum Scissors	85
Injection Cartridge Syringe	86
Tissue Pliers	87
Aspirator	90
Sonde / Irrigating Syringe	91

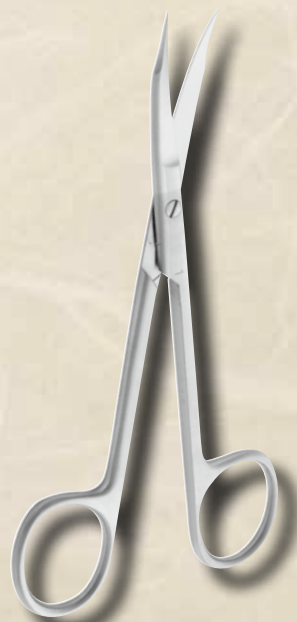


Hemostatic Forceps

Autoclavable up to 135°C

-Mosquito-




Str. ST
No. 29103

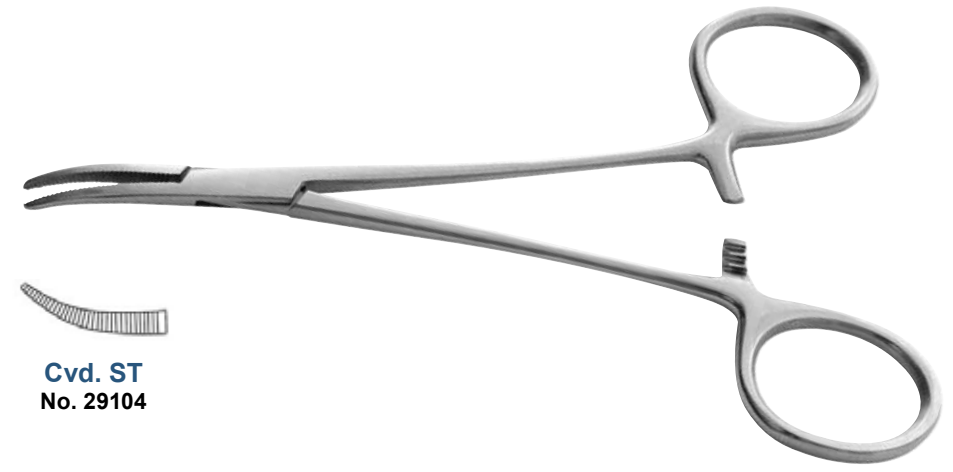


Cvd. ST
No. 29104

Curette

Autoclavable up to 135°C

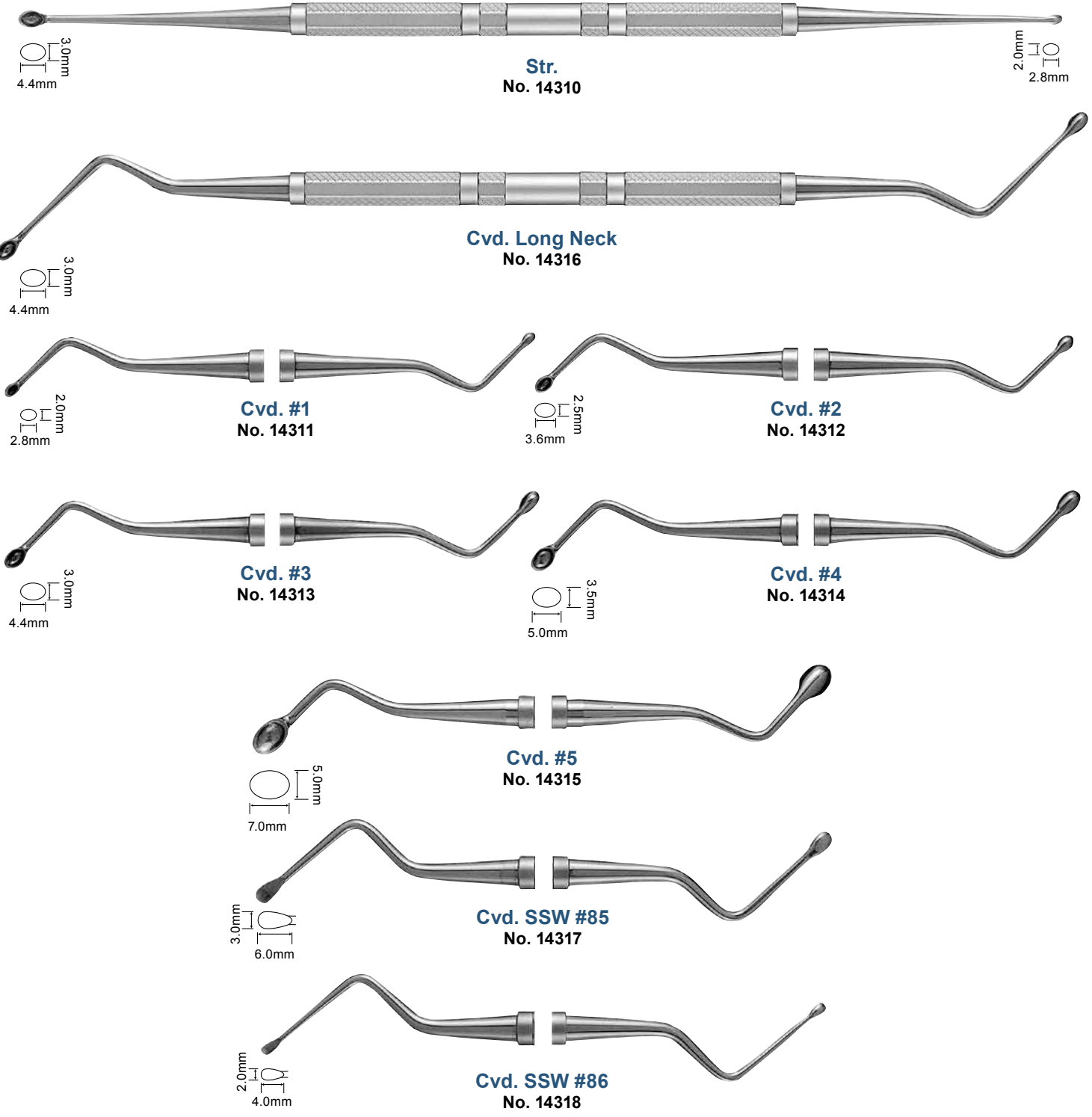

-Windowed Curette-

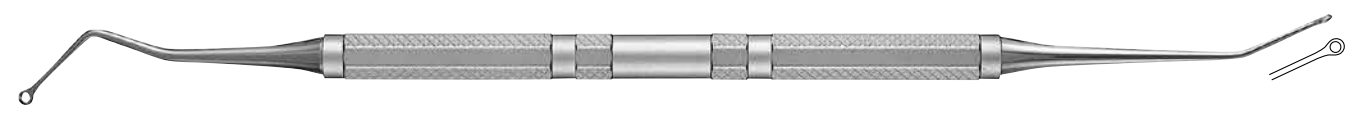

No. 14319

Elevator






-B Handle-



90mm

23mm

Str.

- 
#1S
 Width 4.0mm
 No. 16101
- 
#2S
 Width 3.5mm
 No. 16102
- 
#3S
 Width 3.0mm
 No. 16103
- 
#4S
 Width 2.5mm
 No. 16104
- 
#5S
 Width 2.0mm
 No. 16105

Angled

- 
#20AL
 Width 3.5mm
 No. 16015
- 
#21AL
 Width 3.5mm
 No. 16016



*With an angled shank can insert posterior teeth smoothly.



*It is effective also in the approaching from lingual side of molar part.

Cvd.

- 
#1C
 Width 4.0mm
 No. 16106
- 
#2C
 Width 3.5mm
 No. 16107
- 
#3C
 Width 3.0mm
 No. 16108
- 
#4C
 Width 2.5mm
 No. 16109
- 
#5C
 Width 2.0mm
 No. 16110

for Wisdom Teeth

- 
#1 for WT
 (Slightly Cvd.)
 Width 3.0mm
 No. 16065
- 
#2 for WT
 (Deeply Cvd.)
 Width 3.0mm
 No. 16066
- 
#3 for WT
 (Oppositely Cvd.)
 Width 3.0mm
 No. 16067

*The tip is 20mm longer than regular elevators.
 *To avoid slippery and easy control of blade, they have finger rest.
 *Handle is marked "WT".

Elevator

-D Handle-



29mm

81mm

Str.



#1S Width 4.0mm No. 16121
#2S Width 3.5mm No. 16122
#3S Width 3.0mm No. 16123
#4S Width 2.5mm No. 16124
#5S Width 2.0mm No. 16125

Cvd.



#1C Width 4.0mm No. 16126
#2C Width 3.5mm No. 16127
#3C Width 3.0mm No. 16128
#4C Width 2.5mm No. 16129
#5C Width 2.0mm No. 16130

-D Handle-

NM-Type



#NM1 Width 5.0mm No. 16055
#NM2S Width 4.5mm No. 16059
#NM2 Width 4.5mm No. 16056
#NM3 Width 4.0mm No. 16057
#NM2A (Oppositely Cvd.) Width 4.5mm No. 16058



Root #NM1 Width 2.5mm No. 16554
Root #NM1S Width 2.5mm No. 16555
Root #NM1R Width 2.5mm No. 16556



NM series, with thin and fan-shaped tips In a twisting action, wide edge works effectively.
 *The thin tip can easily get into narrow root space.
 *Its widely designed both ends of the edge effectively work on the root.

[Clinical Movie]



[Clinical Movie]



#NM2 for WT Width 4.5mm No. 16095
#NM2A for WT (Oppositely Cvd.) Width 4.5mm No. 16098

*Handle is marked "RT".

*The tip is 10mm longer than regular NM elevators.
 *To avoid slippery and easy control of blade, they have finger rest.
 *Handle is marked "WT".

Autoclavable
up to 135°C

Autoclavable
up to 135°C

Elevator

Root Tip Pick


-B Handle-

Spear



*With a thin spear-like tip can insert in narrow periodontal membrane cavity easily, and periodontal ligament fibers are cut and it exfoliates.
*Excellent for immediate implantation.

*Since this product is an instrument which cuts a periodontal ligament and carries out extrusion of a tooth by wedge action, the blade is very thin. Therefore, please do not perform operation in which a root is raised by leverage, there are possibility of the breakage of blade and fear of distortion.

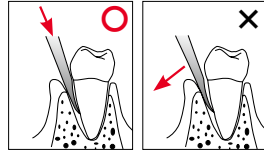

3.7mm



#1S
No. 16141



#1C
No. 16144



#1R
No. 16147



#2
No. 16145




#3
No. 16146

[#1C (Curve)]



*The curve of the tip can approach the "molar" smoothly.

[#1R (Reverse)]

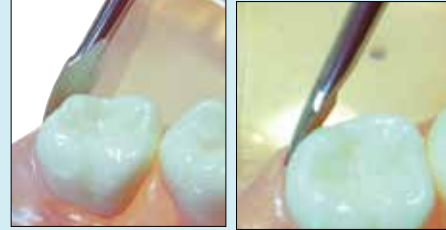


*Inside curve of the tip can approach the end of molar, lingual and palate side.

[#2 #3]



*This curved shape of a shank is suitably effective for "molar" to approach with ease.

*The tip is positioned on the extended line of the handle for handling balance. It makes Dr,s comfortable to keep the balance of operation.

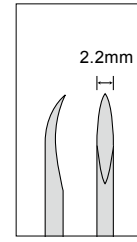

#1S
No. 16085




#2A
No. 16086



#3A
No. 16087

#2A / #3A

*With a angled shank can insert posterior teeth smoothly.

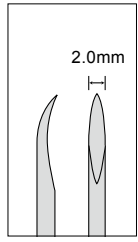


Double End
No. 14601

Root Tip Elevator

*This thin tip of the elevators has easy access to stricture periodontal ligament gaps.

*You can choose one of the most suitable elevators with many types of angles, depending on angle, mouth-opening capacity or condition of remaining roots.

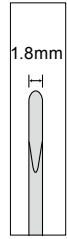

Str.



#1S-0
No. 16810




*The straight type is preferable, which is able to align the line of the tip in that of a front tooth. When it come to molar teeth, it is surely effective for the first broadening of the periodontal ligament gaps as well.

Cvd.




#1S-15
Cvd.15°
No. 16811

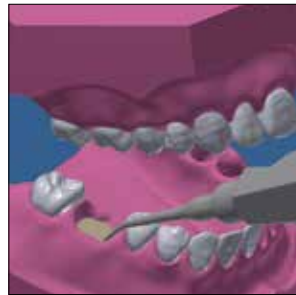

*It is relatively difficult to align the line of the elevator's tip in the line (axis) of a molar tooth. So choose the appropriate one of the 2 curved elevators, depending on the position of the target root and type of patient's mouth.



#1S-30
Cvd.30°
No. 16812

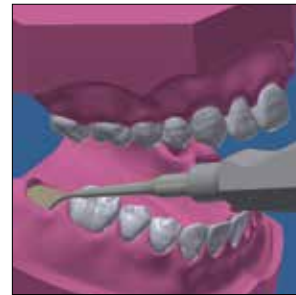

Oppositely Cvd.



#1S-15R
Oppositely Cvd.15°
No. 16813



#1S-30R
Oppositely Cvd.30°
No. 16814



*Choose one of the 2 types of reversed tip elevators, 15 or 30 degree, in case you insert the elevator from the molar, tongue, or roof side.

Angled



#2A
Cvd.30°
No. 16076



*Choose 2A or 3A, especially for those who can't open his mouth widely.



#3A
Cvd.30°
No. 16077



Lower molar L6, distal root



Lower molar L6, mesial root



Lower molar R7



Lower molar L7 mesial root



Periosteal Elevator



Nanaura
No. 14230



#9
No. 15131 *Molt #9* for wide periosteal cutting, retracting or protecting.



#13
No. 15141 *For wide cutting mucossa, cyst or retracting or protecting tissues.



#14
No. 15140



Dissector

*Used to mark on the gum before cut.

No. 30020



**Go your Own way
with your Own choice**



Periosteal Customize System (Not Replaceable Type)

Tip
 *Make your own instrument by combining two tips and handle.
 *Medical treatment becomes more smooth by customized instrument suitable for the technique of the doctor.
 *Price include 2 kinds of any tips.


Handle

-Mucosa (Wide area)-




HA1 HA2 *with Blade HA3 *with Blade HA4 HA5 *with Blade HA6 *with Blade HA7

-Mucosa (Narrow area)-

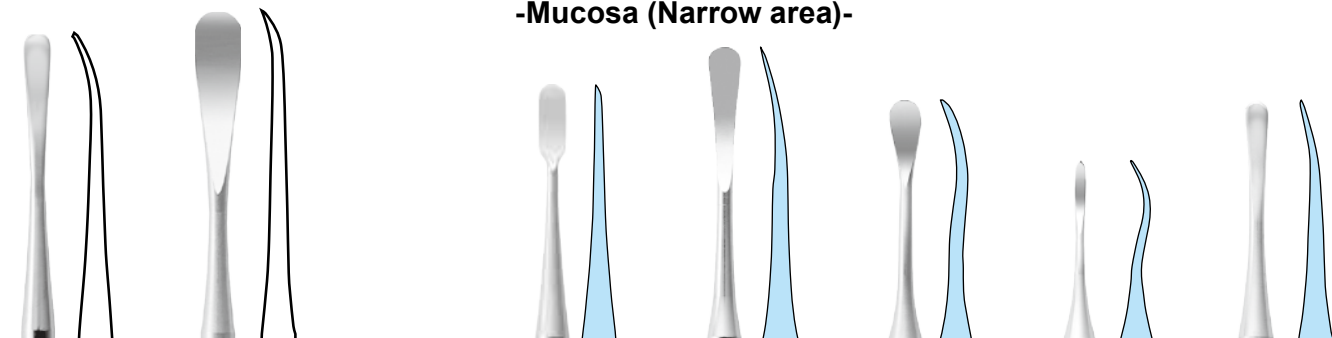


HA8 HA9 *with Blade HB1 *with Blade HB2 *with Blade HB3 *with Blade HB4 HB5 *with Blade

-Interdental Papilla-

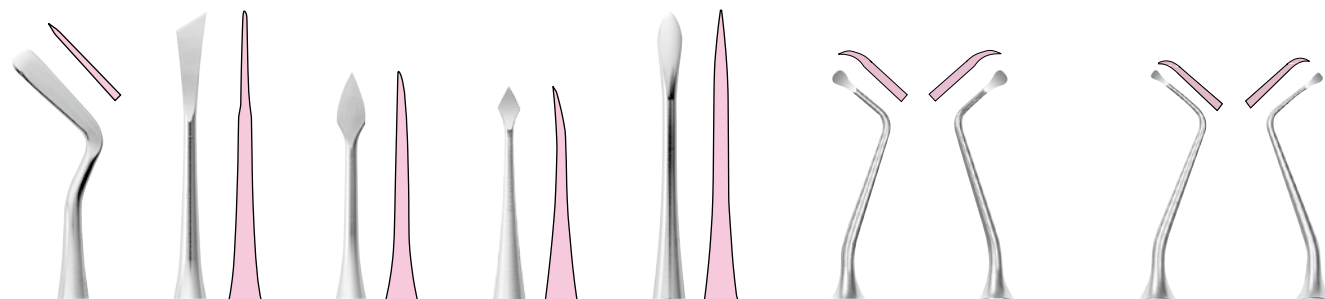

HC1 *with Blade HC2 *with Blade HC3 *with Blade HC4 *with Blade HC5 HC6 *with Blade HC7 *with Blade HC8 *with Blade HC9 *with Blade

-Periosteum-

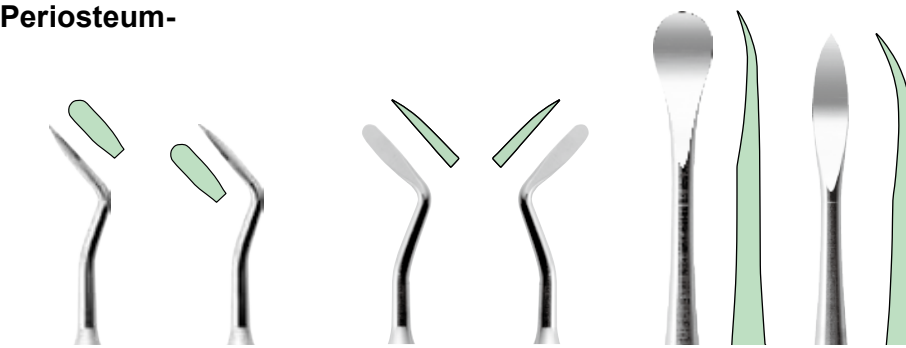

*Backside is flat HE1 *with Blade *Inside is flat HE2 *with Blade HE3 *with Blade HE4 *with Blade HE5 *with Blade HE6 *with Blade

Round No. 15008



*Size / 100mm
 *Handle diameter / ø9.5mm

Round Short No. 15009




*Size / 80mm
 *Handle diameter / ø9.5mm

Octagonal S No. 15003



*Size / 110mm
 *Handle diameter / ø6.5mm

-Gum retractor-

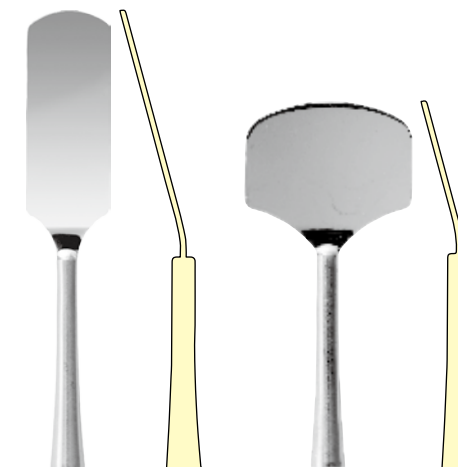

HD1 *Mirror HD2 *Mirror


Next Page Example for Order



Periosteal Customize System (Not Replaceable Type)

Example of Order



*Pattern for gum retractor. The tip of two kinds of size can be properly used depends on operative site.
*It can be used as a mirror.

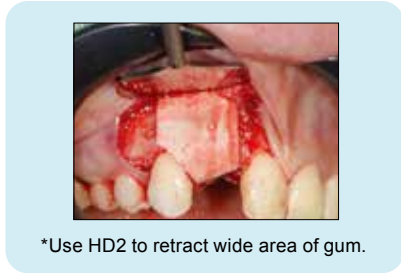


Example of Order

*It is the scalprum by which the angle of four directions was given to the neck.
*They can be respond to all area and ideal for periodontal surgery.

[Clinical Movie]

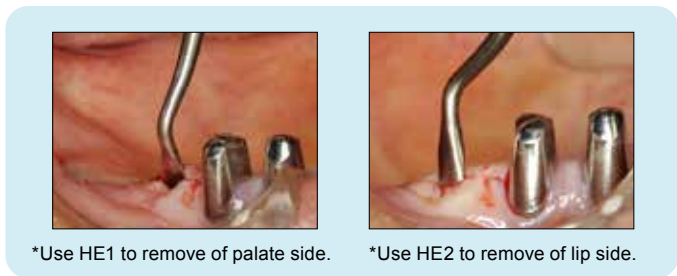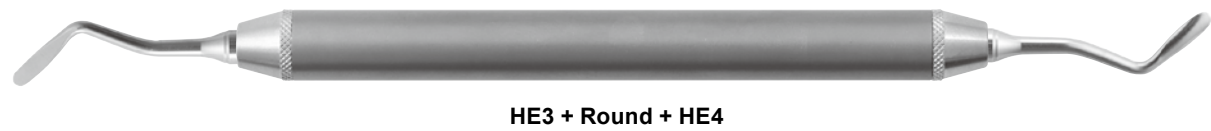

Periodontal Knife

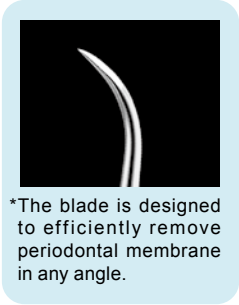

Scalpel Handle





Mallet




Universal
No. 14672

*Weight / 120g
*Material / PTFE (White part)



Hard
No. 15211

*Weight / 220g
*Material / PFA (White part)

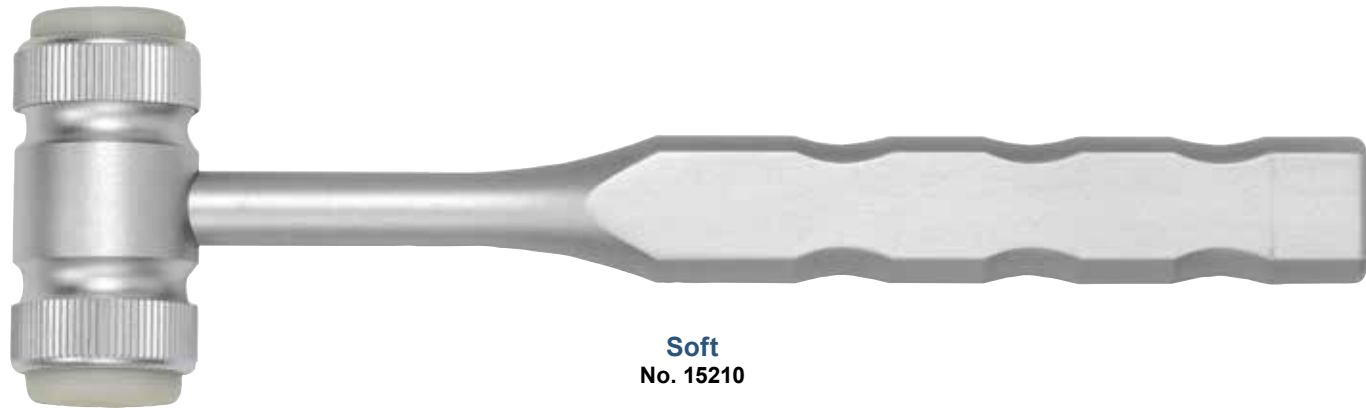

Soft
No. 15210

*Weight / 220g
*Material / Si (Gray part)



Tongue Depressor

No. 14661




Size / 145mm

L-Shaped Retractor



*Suitable for gum valve opening and holding during the surgery operation.

No. 15314

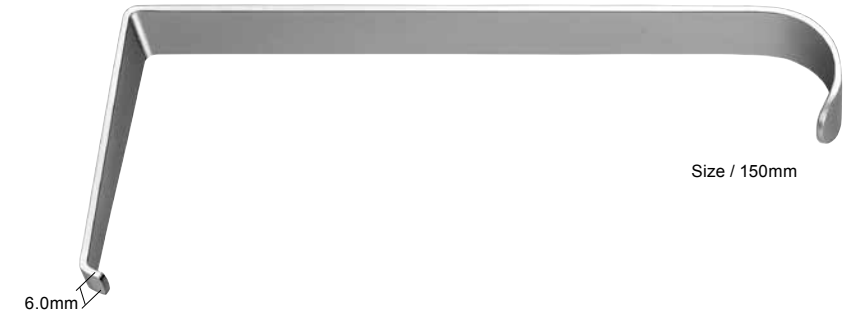

Size / 150mm

6.0mm



Cheek Retractor




No. 14620

Size / 150mm



*Suitable for displace and protect fundal part of periosteal flap.
*Also can be used as Tongue Depressor.



Stemberg L
No. 14657

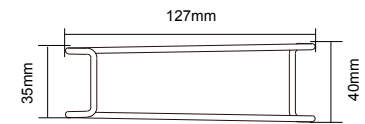

Stemberg S
No. 14658

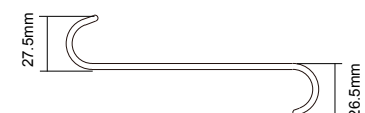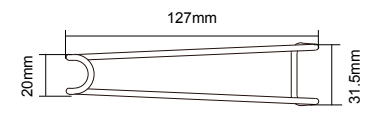




Wound Retractor

-Langenbeck-


*Able to displace efficiently with flat and long hook during surgery operation.
*Even if an assistant holds, a handle is long so that it may not become obstructive and it does not become the hindrance of an operation.



Size / 210mm

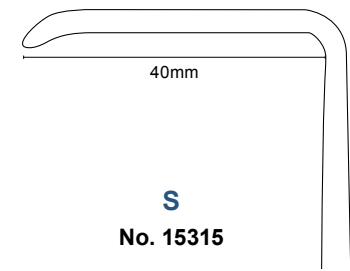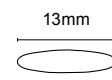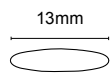

S
No. 15315

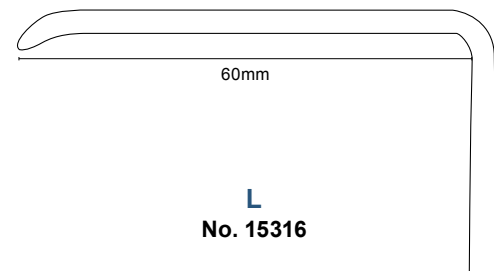

L
No. 15316



Gum Scissors


-Goldman Fox-

#1 ST
No. 29105
*One blade serrated



Size / 130mm

-Iris-

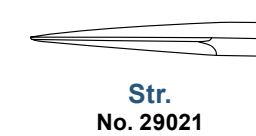

Str.
No. 29021



Cvd.
No. 29022



Size / 114mm



Injection Cartridge Syringe Type-2

-Screw-Type-

**"Light weight" for comfortable operation without getting your hand tired. Approx. 45g ea. (Both No.14101 and 14102)
 **"Suitable length" even for small hand. Max 160mm from top to pushing point. (Both No.14101 and 14102)

[Clinical Movie]



A
Infiltration Anesthesia
No. 14101



B
Conduction Anesthesia
No. 14102



*No.14102 has "spiral screw" to fix cartridge stably.

*Material / PPS (Black part)



**"Rotating handle and push-rod" makes your operation easier from any angle and direction.



**"Resin handle and push-rod" to avoid slipping for stable injection.

Tissue Pliers

-McIndoe-




with Hook
No. 11335

without Hook
No. 11336

-Adson-

*Soft tissue is delicately hold at the time of a suture and an operation.


Str.

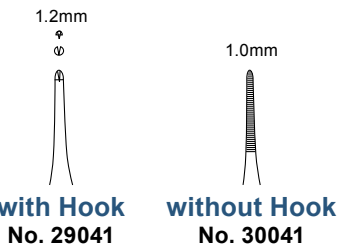

with Hook
No. 29041


without Hook
No. 30041

Cvd.

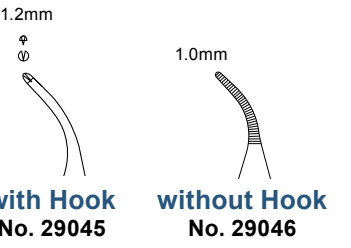

with Hook
No. 29045

without Hook
No. 29046

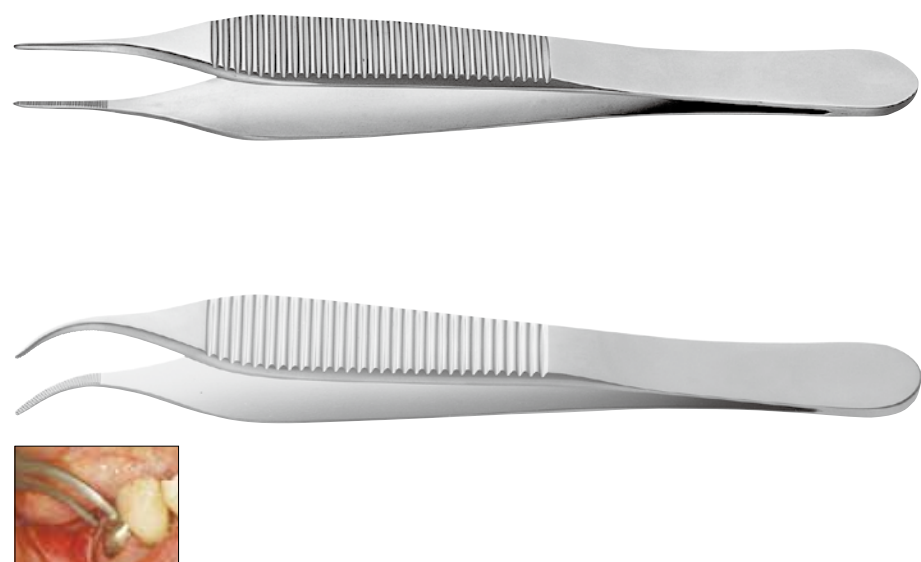


-Graefe-

*Alternate serrated tips will not damage the gum during suturing.

Str.
No. 29043

Cvd.
No. 29044

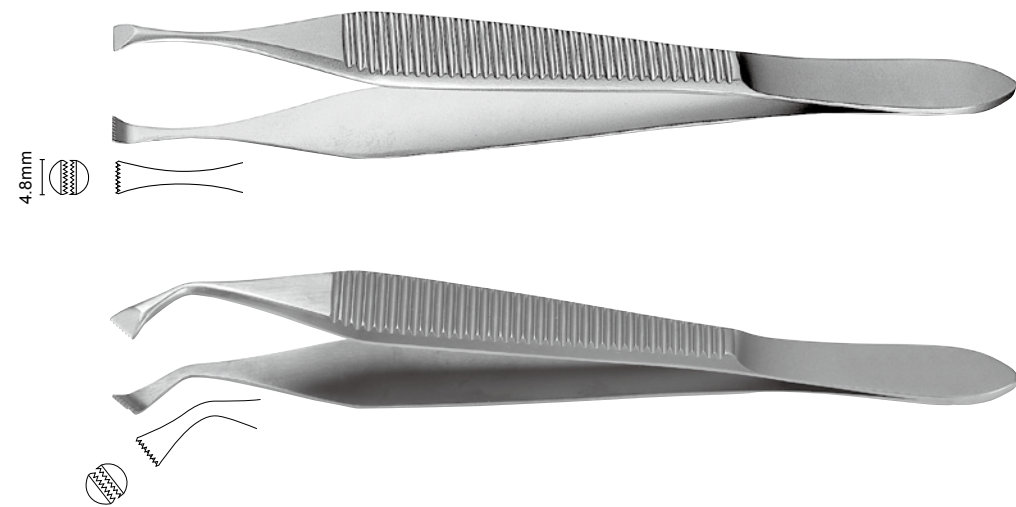

Tissue Pliers

-with Hook-



Str.
No. 11333



Cvd.
No. 11334




#4
No. 15040

-without Hook-




Str.
No. 11331

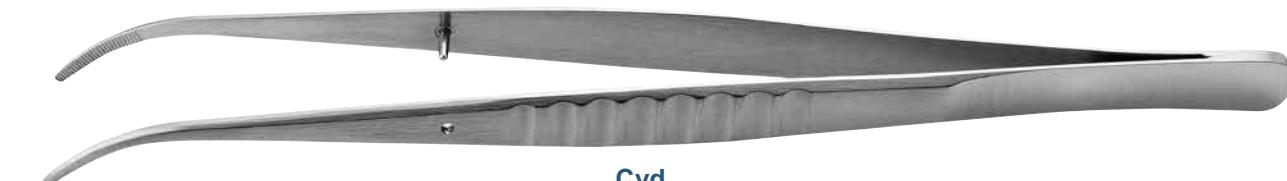

Cvd.
No. 11332

-with Curette-

*With curette shaped tip, able to hold and remove the granulation tissue.

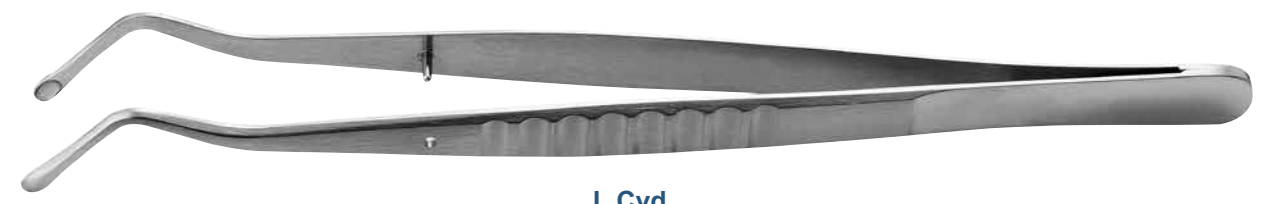

3.1mm
4.5mm

L Cvd.
No. 11324



2.6mm
4.0mm

S Cvd.
No. 11326


-Lucae-

No. 11210




-Suture Pliers-

No. 11216

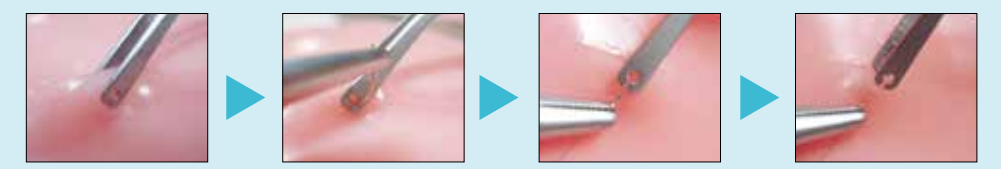


*Prevents gap during suturing by using tip to pinch the gingiva and suture through the hole.

Aspirator

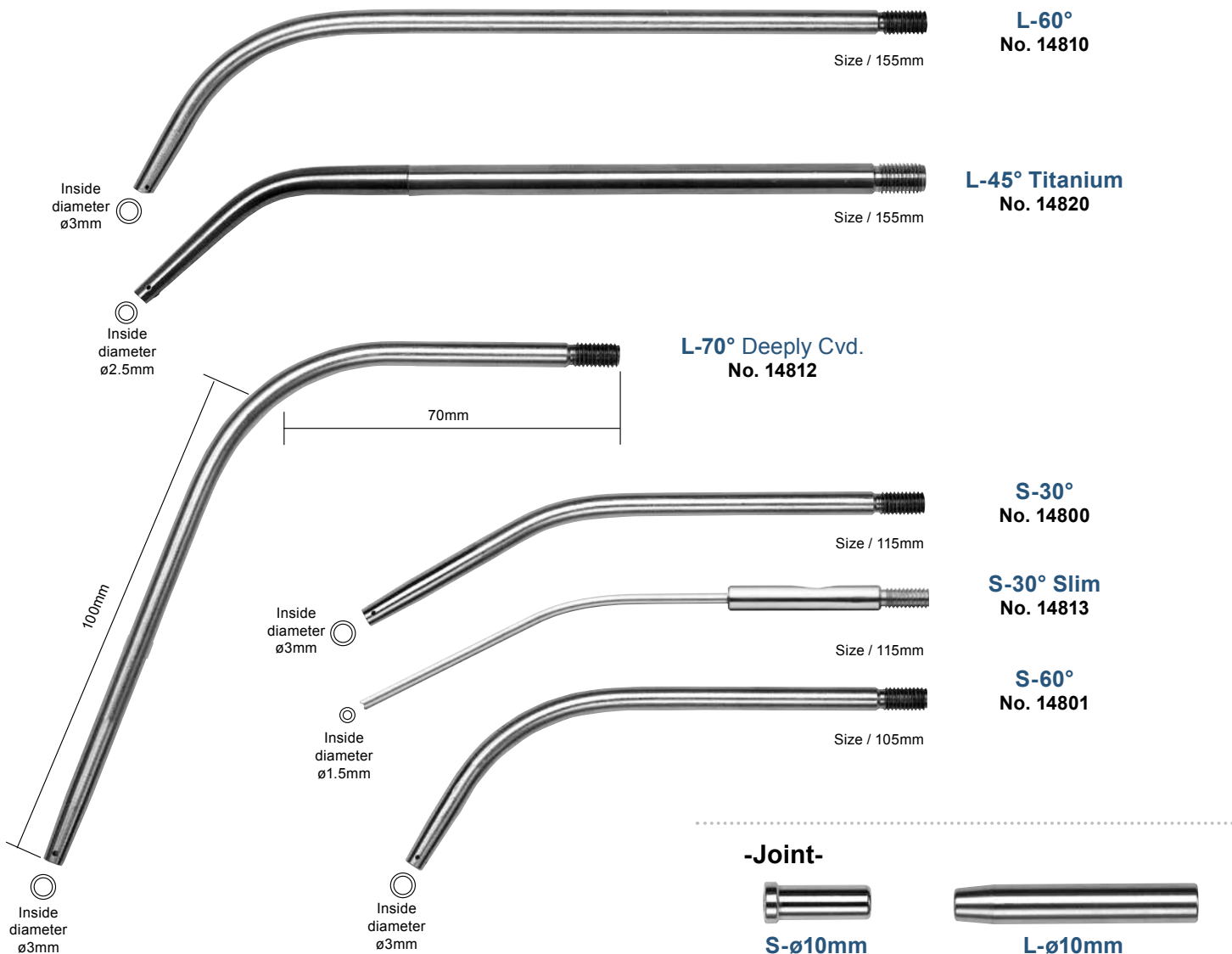

In order to use Aspirator, please purchase main body and joint (1 ea.)

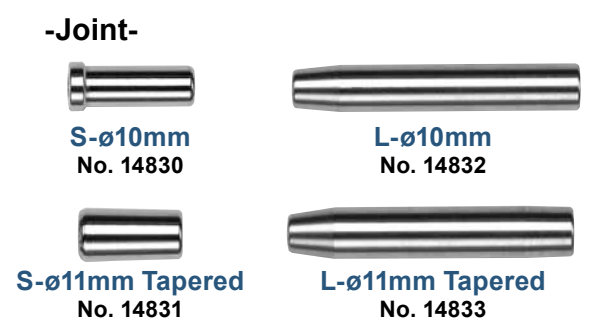

Sold Separately

Adaptor for Saliva Ejector
No. 14840

Autoclavable up to 135°C, Avoid Ultraviolet Sterilization

*Able to connect the aspirator to saliva ejector joint (not evacuation tube).

*Size / 24mm
*Material / Si

Adaptor
No. 13614

Autoclavable up to 135°C, Avoid Ultraviolet Sterilization

*In case the inside diameter of evacuation tube is 16mm, this adaptor is required for use of the aspirator.

Cleaning Brush S (ø 4mm)
No. 13651 (12 ea.)

Not Autoclavable

Size / 200mm

Cleaning Wire
No. 13654 (12 ea.)


Autoclavable up to 135°C

*For "Aspirator S-30° slim", (No. 14813)
*ø1.1mm

Size / 150mm

*Joint ø11 with Tapered is needed for the connection with aspirator.

Sonde

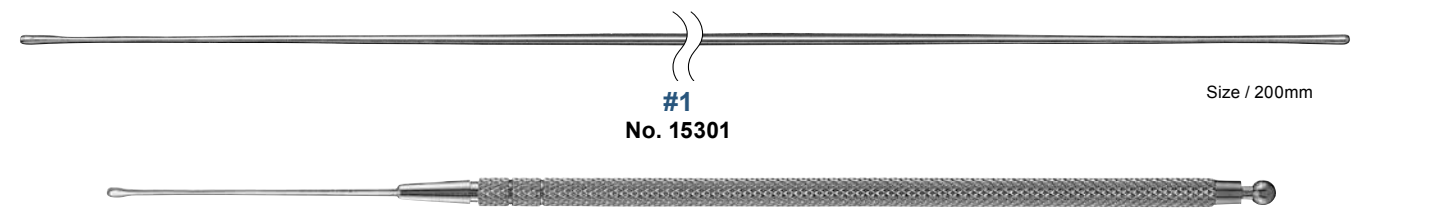


Sonde Holder Set
No. 15020

Contents Holder / Tip

Individual Parts	
No.	Description
15021	Tip
15022	Holder




Irrigating Syringe

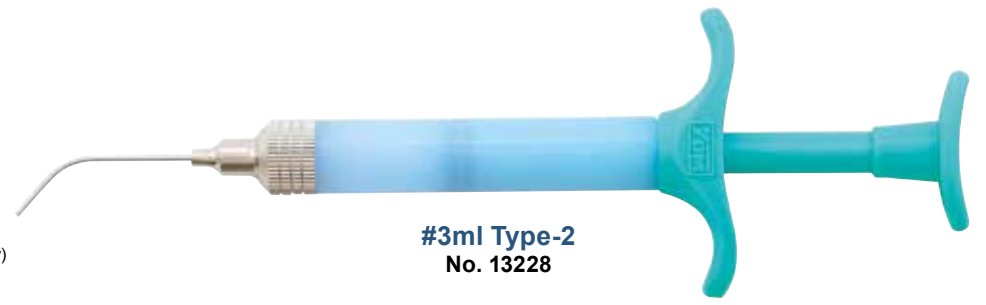

*Material / PFA (Main body)
PP (Trigger)



Jaw Opening Trainer

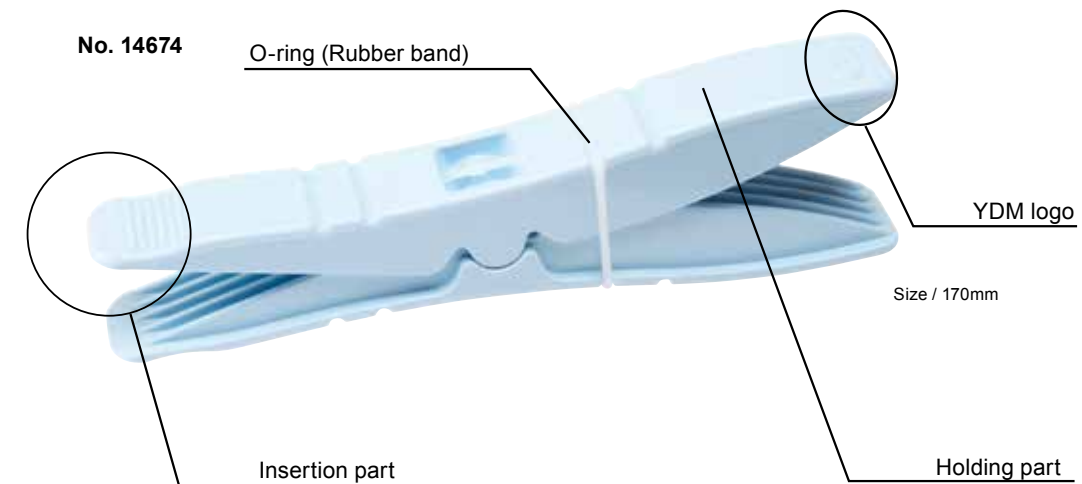

Sold Separately

Silicon Cover (2 ea.)
No. 90424

Autoclavable up to 135°C, Avoid Ultraviolet Sterilization

*Bite pressure / up to 490N
*Material / PC, ABS (Main body)
Si (O-ring)

*Structure automatically opened by the power of the O-ring hung on the opposite side.

-How to Use-

*Photo by: Tokyo Medical Dental University (Temporomandibular Joint Clinic)

- 1) Close an insertion part (the side which does not have a YDM logo) by hands and inserts into a mouth.
- 2) Adjust the dimpled area to the canine.
- 3) Grasp at the holding part lightly and an insertion part will open by the power of O-ring.

