



Bone Grafting Product	94
Sinus Lift Instrument	96
Periotome / Elevator	97
Measuring Device	98
Titanium Instrument	98
Resin Color Probe	99





Bone Scraper Type-2



Size / 155mm



*Clear body allows how much you scraped own born inside to be visible. *Suitable angle of the tip part allows easy to scrape bone.

*Collect own born like shaving allows less surgical invasive.

*Capacity / approx. 1cc *Width of blade / 6.4mm *Material / TPX

Bone Mill



-ex-2-



Graft Carrier

*Used for transfering of a autologous bone or an artificial bone under implant treatment.



No.	Description	Diameter (mm)	Capacity (cc)
30462	ø2.9	ø2.9	0.02
30463	ø3.4	ø3.4	0.03

Dissector



*Used to mark on the gum before cut.

No. 30020







Sinus Lift Instrument



- *Smooth edges to prevent Schneider's membrane from damages.
- *Easy control angle design.
- ${}^{\star}\text{Suitable}$ for searching bone deeper inside maxillary sinus.

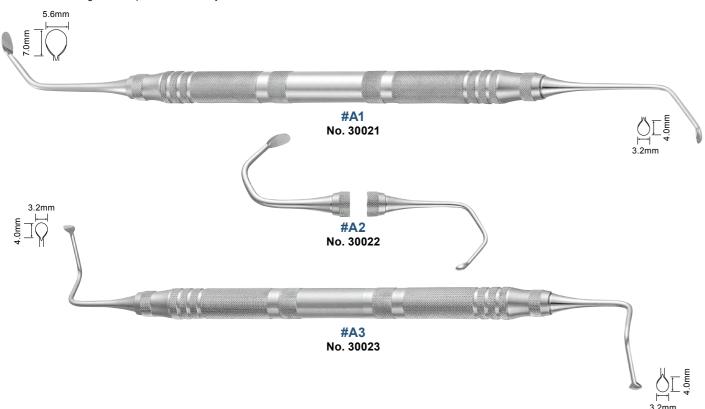

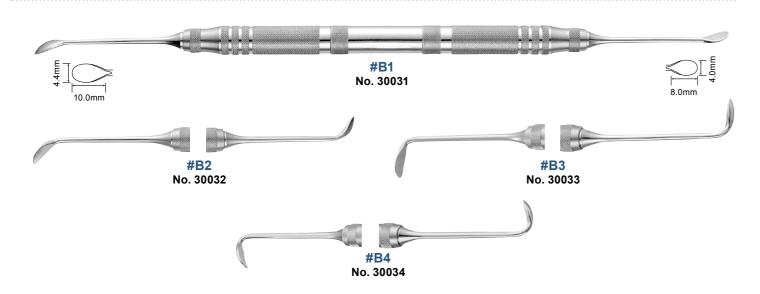



Periotome



- *Able to operate extraction with minimum damage to surrounding periodontal tissues.
- *Excellent for immediate implantation.





-Handle-



P-A (ø8mm)


30016	Point #1
30017	Point #2
30018	Point #3
10201	P-A Handle

Elevator






*With a thin spear-like tip can insert in narrow periodontal membrane cavity easily, and periodontal ligament fibers are cut and it exfoliates.
*Excellent for immediate implantation.

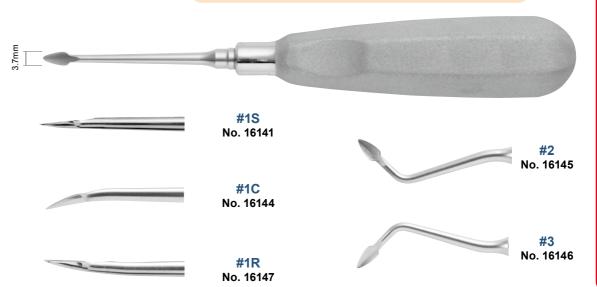

*Since this product is an instrument which cuts a periodontal ligament and carries out extrusion of a tooth by wedge action, the blade is very thin. Therefore, please do not perform operation in which a root is raised by leverage, there are possibility of the breakage of blade and fear of distortion.

Description

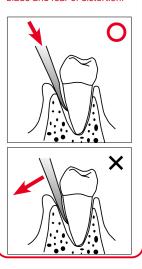



*The curve of the tip can approach the "molar" smoothly.



*Inside curve of the tip can approach the end of molar, lingual and palate side.



- *This curved shape of a shank is suitably effective for "molar" to approach with ease.
- *The tip is positioned on the extended line of the handle for handling balance. It makes Dr,s comfortable to keep the balance of operation.



Parallel Meter







X-Ray Gauge



Contents Sleeve (2 ea.) / Gauge (2 ea.)

No. 21157



Sleeve



Gauge (1mm by 1mm with hold line every 5mm)

*Suitable for root canal treatment, accurate examination before placing an implant, checking of calculus quantity and its area.



W33 x D45mm (Gauge) *Material / Si (Sleeve)

Periodontal Probe

*Size / W38 x D48 x H8.5mm (Sleeve)



*Titanium




Periodontal Color Probe Tip (Replaceable Type)



-Tip (5 ea.)-

- *A probe made of the resin prevents implant and a prosthesis being damaged.
 *Black and yellow color combination offers easy distinction of the scale.

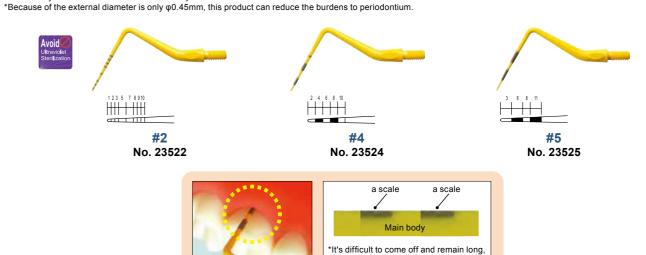

*Material / PPSU





because of the laser marking permeates

inside the resin body.



*Material / POM

99

Mouth Prop








*Able to keep the mouth open stably even for surgery which takes some hours especially under sedation.


-with Tongue Guard-


Contents Mouth Prop / Tongue Guard Right, Left



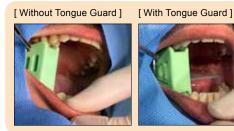



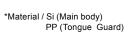

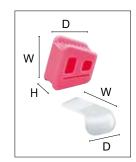

Description		Size (mm)		
		W	D	Н
Mouth Prop	L	40	35	17
	S	35	29	16
	SS	30	25	16
Tongue Guard	L/S	52	34	-
	SS	39	24	-

Almighty Mouth Gag



*Able to keep the mouth open stably even for surgery which takes some hours especially under sedation.

Contents Main body, Tip (4 ea.)





^{*}Bite pressure / up to 1400N *Material / Si (Tip)