




184 185 192
192
107
197
200
203
207
207
208
209

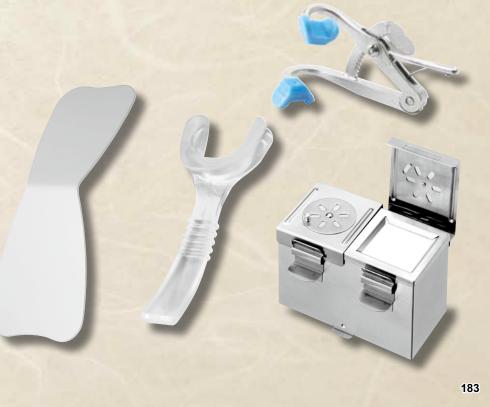

Mouth Prop

*Able to keep the mouth open stably even for surgery which takes some hours especially under sedation.

-with Tongue Guard-

Contents Mouth Prop / Tongue Guard Right, Left






Description


Mouth Prop

Tongue Guard


utoclavab up to 135°C ()

Cheek Retractor

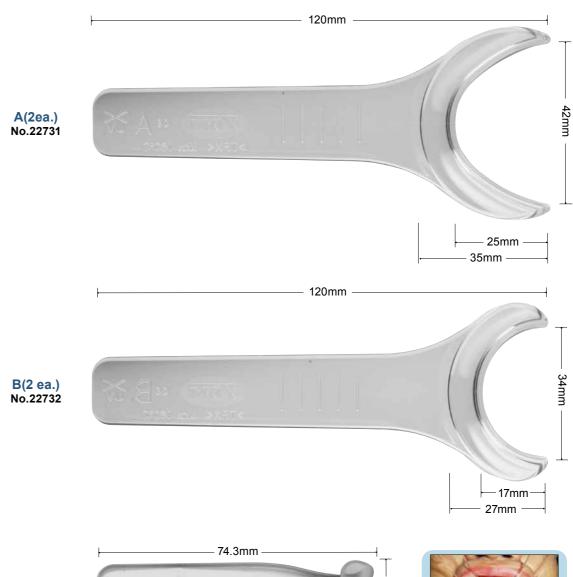

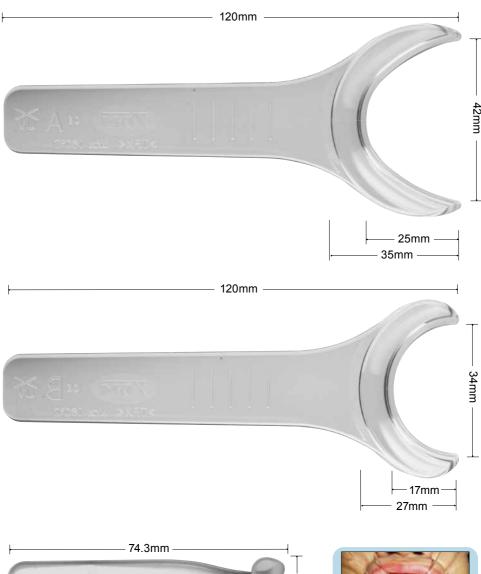

15mm

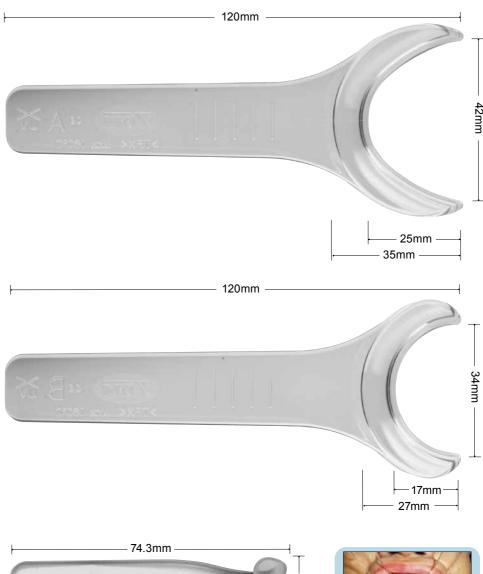

*Material / Si (Main body) PP (Tongue Guard)

Almighty Mouth Gag

*Able to keep the mouth open stably even for surgery which takes some hours especially under sedation

SS

No.14678

Contents Main body, Tip (4 ea.)



D

No.14670

D

W



Size (mm)

D

35

29

25

34

24

н

17

16

16

W

40

35

30

52

39

L

S

SS


L/S

SS

C(2 ea.) No.22733



D(2 ea.) No.22734

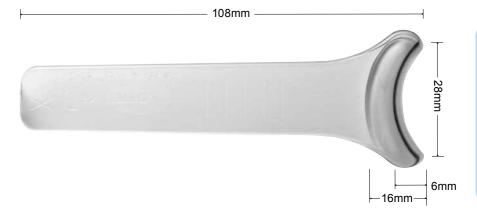

*Material / TPX


*Bite pressure / up to 1400N *Material / Si (Tip)





photography, possibly used as lip retractor. With photographic mirror occlusal.

Cheek Retractor

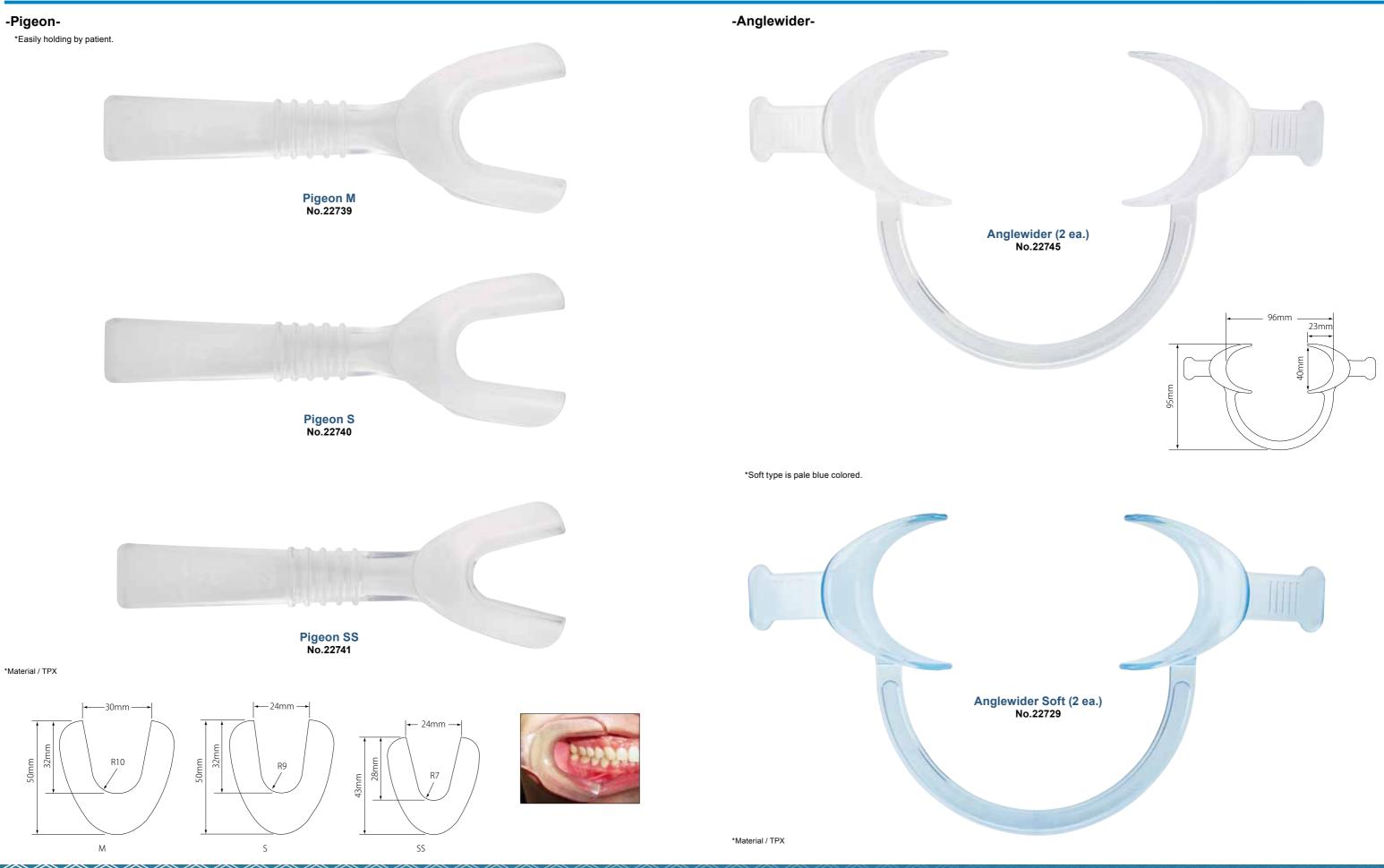



Cheek Retractor

-Anglewider-


*Easy insert by 1 hand with soft material keeping enough mouth-opening. *Able to keep the clear view without any assistance. *Soft material for patient's comfort. *Able to expand the view by pulling its wings.




Cheek Retractor

-Anglewider-

*It opens corners of the mouth, and allows clear views without auxiliary personnel. *The inner wings retract buccal mucosa, and thereby oral vestibule is clear.

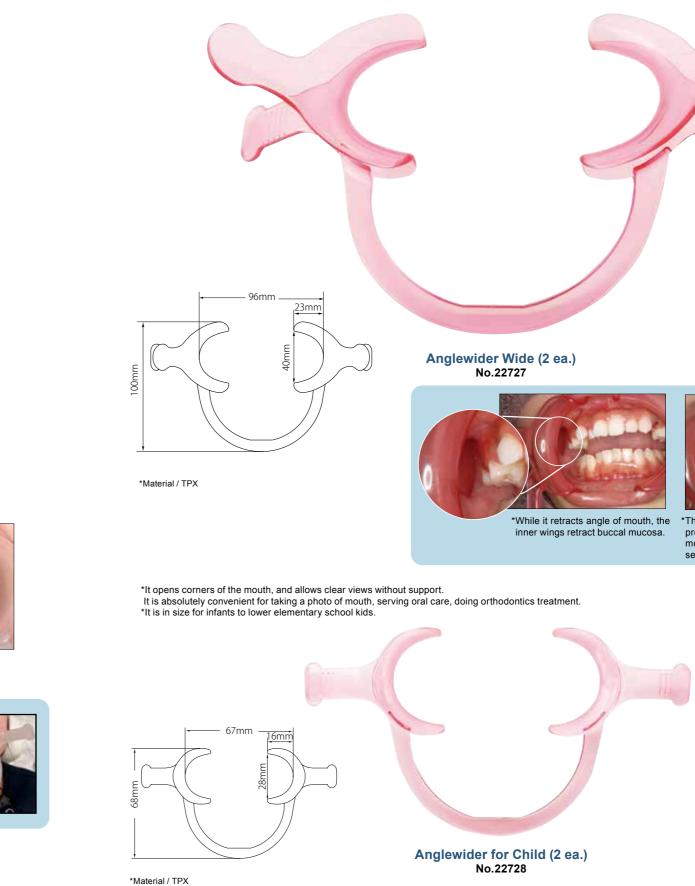

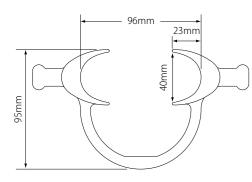

*Material / PP : Deflection temperature 120°C

Anglewider One Hand (2 ea.) No.22726

















*The wings retract oral vestibule, which provides you with comfortable environment to bond orthodontics brackets, serve oral care and fill resin as well.

Cheek Retractor

-Doublewider- *Retract cheeks, lips, and tongue for maximum visibility and patient's comfort. *Suitable for dental aesthetic instruction, inoffice bleaching and bracket wearing in orthodontic treatment.

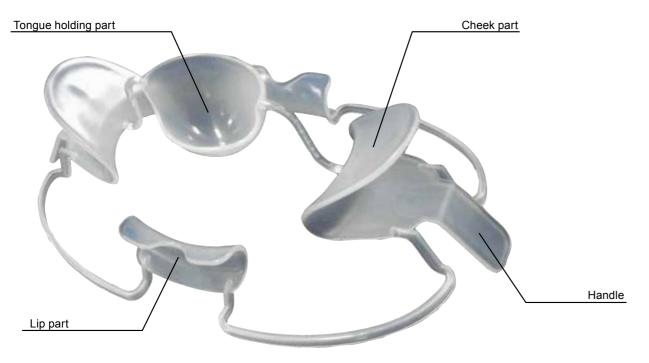

*Material / PP : Deflection temperature 120°C

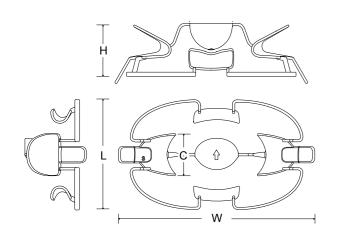

No.	Description		Size	(mm)	
INU.	Description	W	L	Н	С
22750	Doublewider L (2 ea.)	171	96	54.5	38
22749	Doublewider S (2 ea.)	161	90	50.5	35
22725	Doublewider SS (2 ea.)	157	80	42	25



*Tongue holding part can be removed if you don't need it.



 Make sure the up and down by the arrow on the tongue holding part, twist the both handle parts and insert the right and left cheek parts into the mouth.










*This dental-arch shape retract the buccal mucosa of molar part.

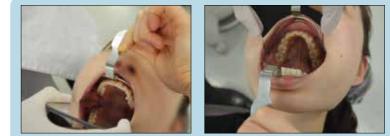

*Patient can hold the handle while operation and the handle is designed not to bother photographing.
















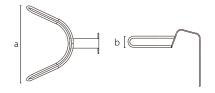

No.	Description	Size (mm)	
INU.	Description	а	b
22723	Deep Retractor SU Adult	80	12
22724	Deep Retractor SU Child	67	9




Photographic Mirror

-Stainless Steel-

*Use of stainless steel avoids breakage, allowing long term use of the mirror. *Sterilization pouch included.

*In sterilization,put mirror into pouch without disposal sterilization pack.

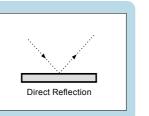

Rodium Coating (With "Rodium" marking on the back side) (1) *Reflection rate / About 85%

*Without Rodium Coating (2) *Economy Type *Reflection rate / About 65%



60mm

*Bended design for keeping fingers out of your clinical photographs.



70mm

60mm






*"Occlusal Wide LL" and "Occlusal LL" Mirror is ideal for taking a photpgraph of "Margin" and "Buccal mucosa"even the patient has wide dental arch.

Buccal No.13573

ST Buccal No.13263

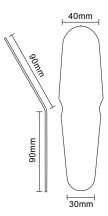




35mm



Photographic Mirror



Comparison of YDM Photographic Mirrors

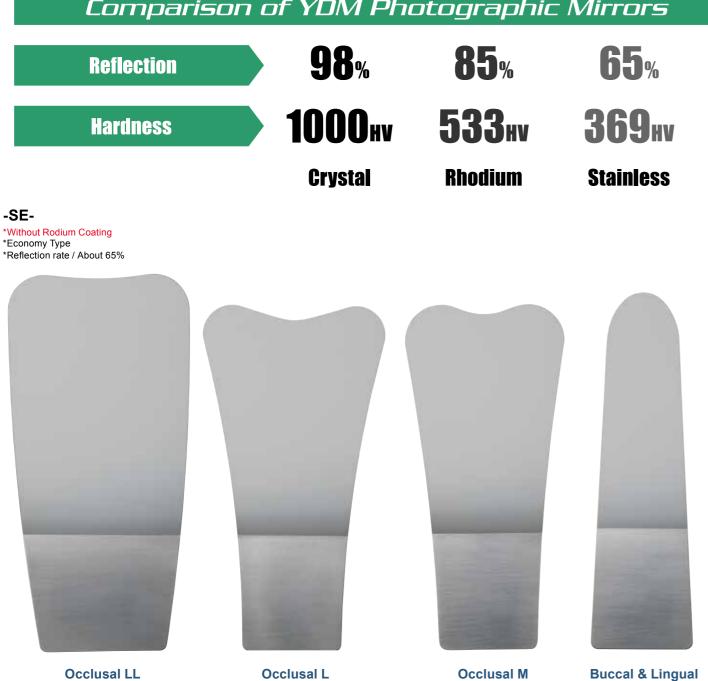


Occlusal LL No. 13274


155mm

75mm

Occlusal L No. 13272

74mm

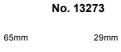

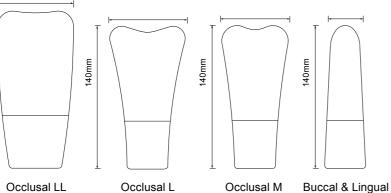



No. 13271

-Made of Japanese glass-*Dielectric multilayer coating

Photographic Mirror Crystal

- *Reflection rate / About 98%
- *Light weight *One sided




Occlusal L


*Thickness / 3.0mm

14 0





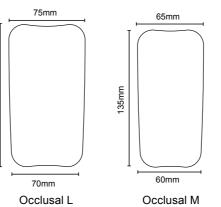

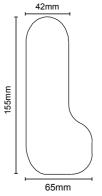





Occlusal M No. 13283

Buccal & Lingual No. 13281



その 他 診療 参 Other Chair-side Products

Contraster

*Teflon coated aluminum.


*Suitable for checking before and after treatment, like filling composite. *Designed to fit the shape of dental arch.









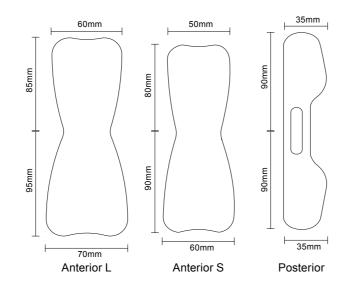

Bur Guard


*Color-coded design allows organization of burs according to their type, treatment method, length, etc. *The cover will not allow burs to fall out.

*Also used for storage for engine reamers.

No.	Description
13471	Yellow
13472	Orange
13473	Pink
13474	Blue
13475	Green



*Storage Capacity / 12 burs. *Holding length / from 17 to 37mm (When using the long cover / from 32 to 52 mm) *Size / W69 x D10 x H41.5 to 56.5mm *Material / POM *Sterilizable with liquid sterilizing solutions.

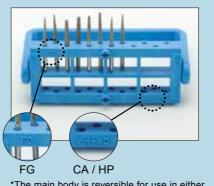



*The main body is reversible for use in either FG burs or CA / HP burs.

*The height can be adjusted to three levels.







*Each "Bur Guard" can be connected.

Sold Separately




No.	Description
13481	Yellow
13482	Orange
13483	Pink
13484	Blue
13485	Green



*Also able to hold laboratory hand piece burs when using the long cover.

Bur Stand

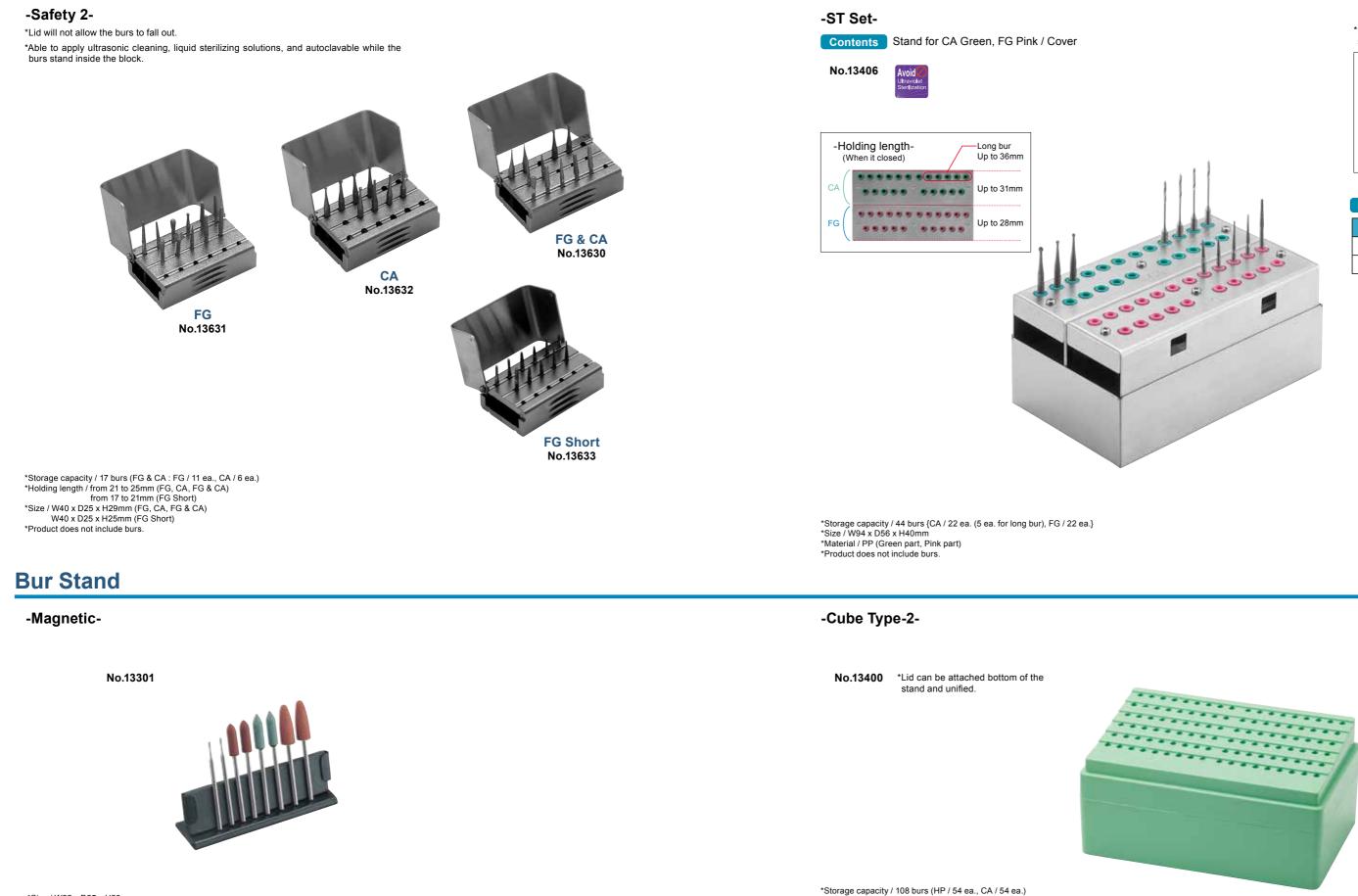

*Holding length (When it closed) / up to 48mm (HP) up to 27mm (CA) *Size / W100 x D63.5 x H56.5mm

*Size / W82 x D25 x H28mm *Some burs will not be hold by magnet. *Product does not include burs.



*Lid can be attached bottom of the stand and unified.



Individual Parts

No.	Description
13407	Stand for CA Green
13408	Stand for FG Pink






One-Touch Cotton Box



Cotton Dispenser

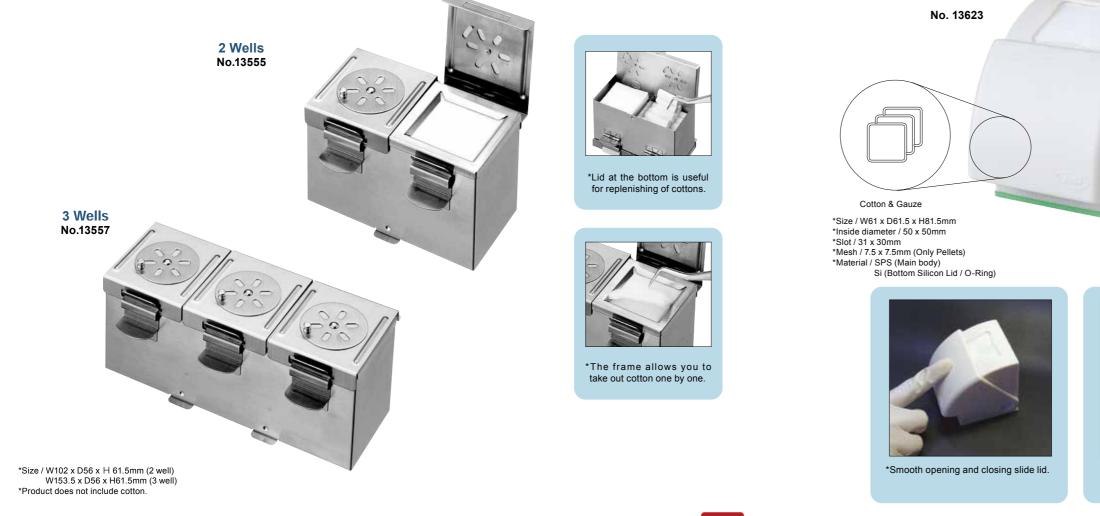

One-Touch Cotton Box Type-2

-for Alcohol-

S No.13621

L No.13622



*Size / ø80 x H57mm (S) ø103 x H82mm (L)

utoclavable ^{ip to} 135℃


Phial Stand One-Touch

*Easy opening of lid by just flipping up the cover.

No.13520




*Size / W100 x D52 x H55mm *Material / Si (Stand base) *Dimension for medicine bottle space. *Product does not include medicine bottle.

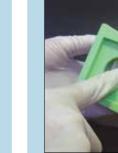


*If cottons decrease, push up the bottom and maintain the position to take out easily.



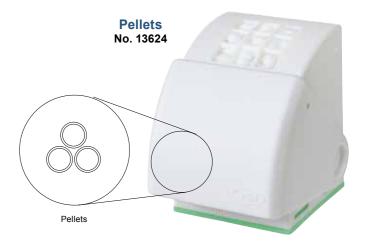





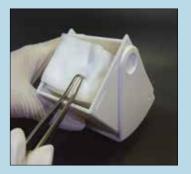

*Since it supplies from a bottom, it can supply without taking out the remaining cotton



Void Iltraviolet



Z-1 eco Fine Liquid

No.13538



High detergency *Alkaline ph works exelently well against protein and fat content. This liquid is efficient for any metal materials due to the special blended component.

Stable detergency *Z1 doesn't matter the temparature of dilution water.(enzyme is not included)

Protect surface * Colloidal particle wraps removed dirt and prevent reattachment.

Rust prevention * Z1 has excellent rust prevention power thanks to its sequestering agent.

Good cost performance and compact bottle *200 fold dilution. Purified water is recommended in dillution. *Easy to meaure by one hand.

Less damage on environment *Z1 is refined from rosin oil , plant based detergent,. *Biodegradation rate is as high as 97.5% ,Z1 can be drained after cleaning

[Form and Component] ■Form : 500ml concentrate(will be 100L by 200 fold dilution) Component : Alkaline surfactant as plant based detergent with rust prevention material *Anionic surfactant, Nonionic surfectant *norganic builder(silicate, carbonate) *Organic builder(carboxylic acid etc..)

*Please follow and confirm the relevant laws or regulations of your country for import.



Sterilization Cassette

-Type-A-

Contents Main body / Single End Clip (3 ea.)

*Stackable cassette. *Clips are used for holding / fixing single-end instrument.

No.24004



*Storage capacity / 8 ea. *Size / W185 x D120 x H25mm *Material / PPS

-Type-B-

Contents Main body / Single End Clip (3 ea.)


*Stackable cassette *Clips are used for holding / fixing single-end instrument.

No.24009

Color Ring (White)



*Storage capacity / 5 ea. *Size / W185 x D80 x H25mm *Material / PPS

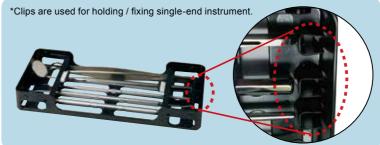







Sterilization Case







Bur Pot Set

Contents Pot / Basket (2 ea.)

No.13304





utoclavable p to **135°C**

Autoclavable ^{up to} 135°C

General Tweezers






*Size / ø75 x H55mm (Pot) ø70 x H45mm (Basket)

Bur Basket



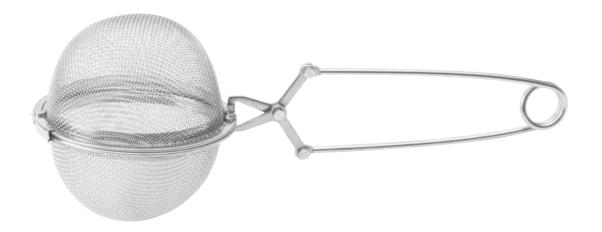

*Size / Inside diameter ø90 x H135mm *Product does not include instruments.

Instruments Stand

No.13547



*Size / ø60 x H66mm



23cm No.13414 Size / 230mm





Color Tape

*Tape is protected by plastic case with open and close lid for sanitary purpose.

No.	Description
13661	Red
13662	Yellow
13663	Blue
13664	Green
13665	White
13666	Black





*Material / Si

^{ip to}**121°**C

Cleaning Brush

Sterilization Cap

*There are two kinds of size.

*Protect cutting edge during sterilization, disinfection and storage.



*Material / Si

*1 roll / about 3m


*Width / 3mm *Size / W8 x D72 x H44mm

*Material / Polyester : Deflection temperature 125°C



No. Description 13523 Red (20 ea.) 13524 Green (20 ea.) 13526 Blue (20 ea.) 13527 Black (20 ea.) 13528 Yellow (20 ea.) 13530 White (20 ea.) 13671 Pink (20 ea.) 13675 Violet (20 ea.)

Outside diameter / ø10mm -L-Inside diameter / ø8mm Size / 5mm

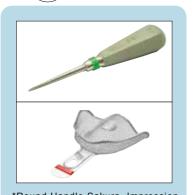

No.

13532

13534

13533

13682

13531

13684

13535

13688

Description

Red (10 ea.)

Green (10 ea.)

Blue (10 ea.)

Black (10 ea.)

Yellow (10 ea.)

White (10 ea.)

Pink (10 ea.)

Violet (10 ea.)

*Round Handle Sakura, Impression Tray, Elevator etc.

Cleaning Wire

ʻø1.1mm

No.13654 (12 ea.)

*For "Aspirator S-30° Slim".

S (ø4mm) No.13651 (12 ea.) *******


L (ø12mm)

No.13653 (12 ea.)







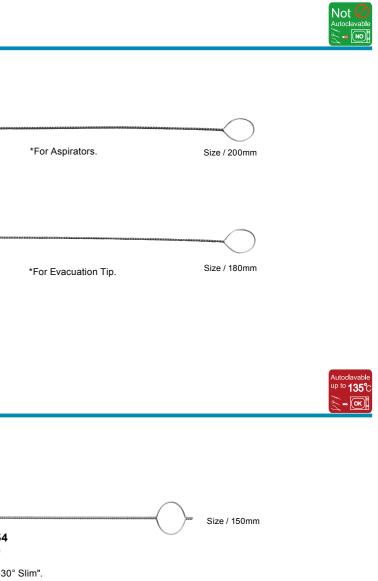



Instrument Oil



*For rust prevention and sharpening.

- $\ensuremath{^*\text{Please}}$ put a liitle of this oil before autoclave for rust prevention.
- $^{\ast}\mbox{This}$ oil works excellent for the joint part of pliers and forceps as well.
- *This oil also can be used for sharpening scaler and curette.
- *Please follow and confirm the relevant laws or regulations of your country for import.

No.13093



20mL