

Extracting Forceps Claw




Patent Pending

4 claws perform excellent hold and dynamic balance. Moreover, the serration of inside the tip makes holding-power tight and strong.



<Evolved Extracting Forceps>

3 or 4 claws always hold any tooth cervix dynamically well. Therefore this forceps are suitable for any shape of tooth cervix.

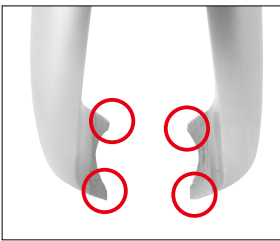

<Rough Surface inside the Tip>

The fulcrums of claws and the edge of rough surface inside the claws make holding power more stronger. Therefore this forceps prevents Dr. from putting too much power on this forceps and this efficiency makes it possible to extract such as fragile tooth or damaged one.

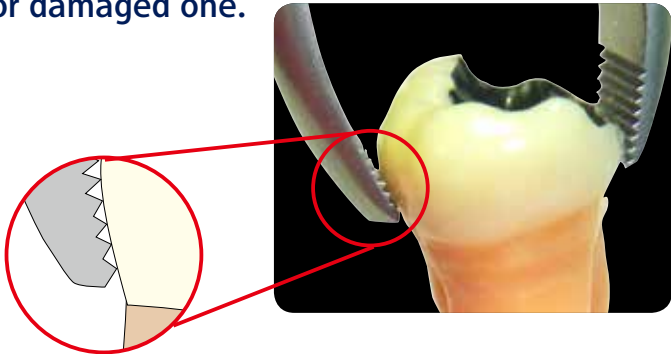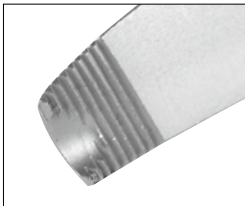


	Standard / used for			
Claw	#1	#4	#10S	#27
	Anterior	Premolar	Upper molar	Lower molar
Claw for milk tooth	#An	#Po-U		#Po-L
	Anterior	Upper posterior		Lower posterior

<Fine Adjustment>

by spring

Thanks to this spring, the range of movement is very stable during an operation. It also means "easy control". Dr. doesn't need to put his/her finger in the space of between handles for adjusting the range of movement like conventional extraction forceps.

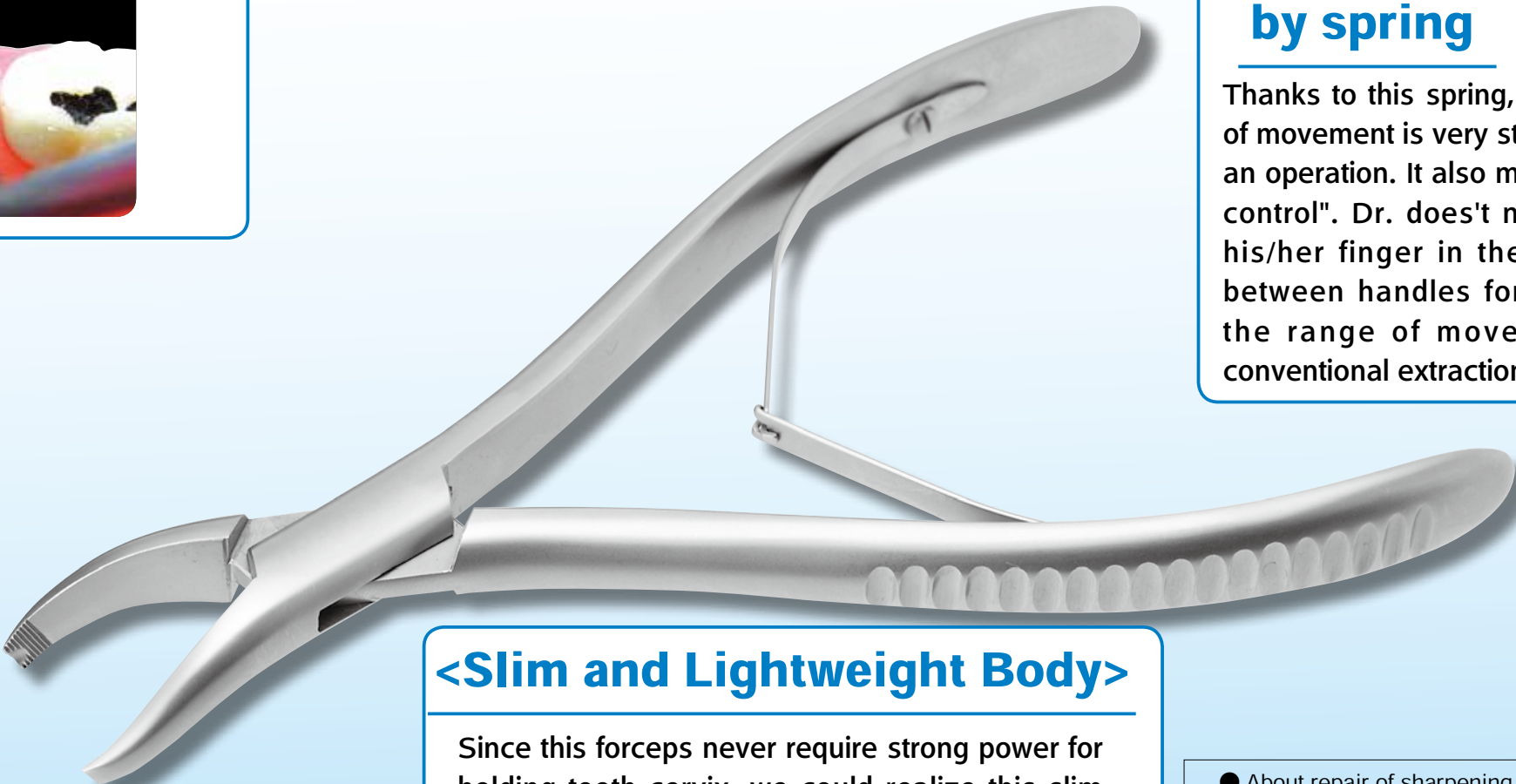


<Slim and Lightweight Body>

Since this forceps never require strong power for holding tooth cervix, we could realize this slim and lightweight body. Excellent comfortable for both men and women Drs.

● About repair of sharpening the edge
We could sharpen the edge which has been got worn.
Please ask the dealer you purchased.

Extracting Forceps Claw

Pediatric Extracting Forceps Claw

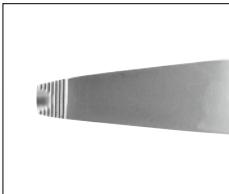


[Actual Size]




#1
No.17581

Size / 162 mm

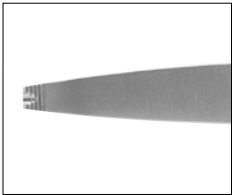

[Actual Size]

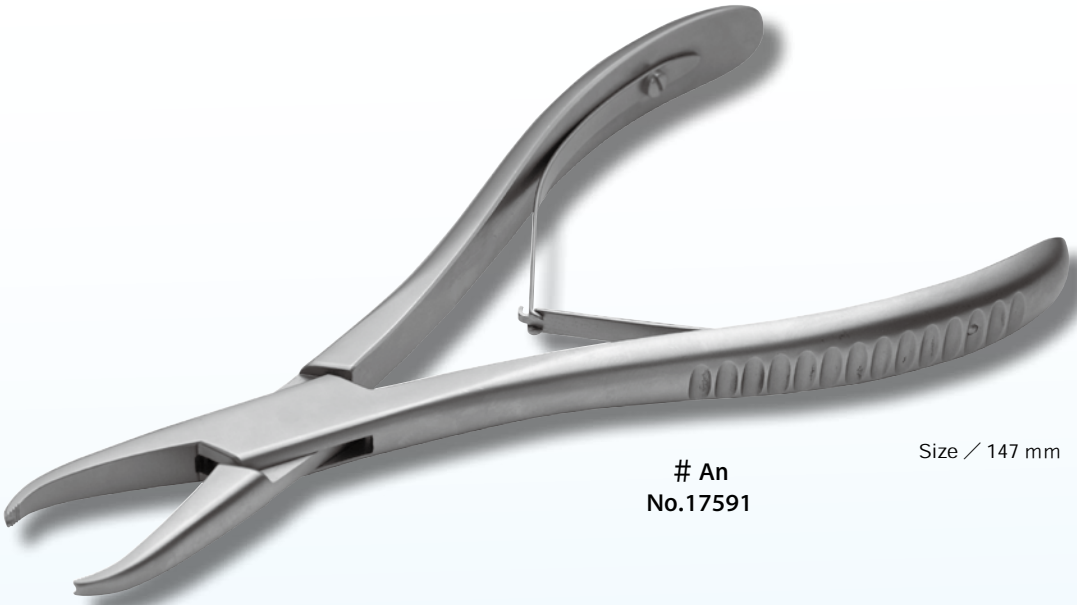


An
No.17591

Size / 147 mm



[Actual Size]

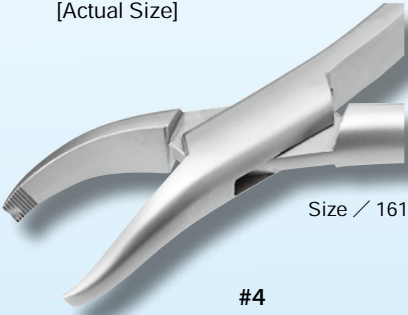


#4
No.17582

Size / 161 mm



[Actual Size]

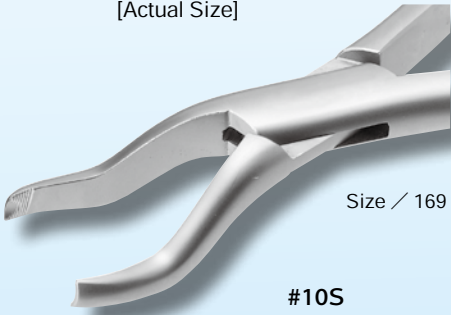


#10S
No.17583

Size / 169 mm



[Actual Size]

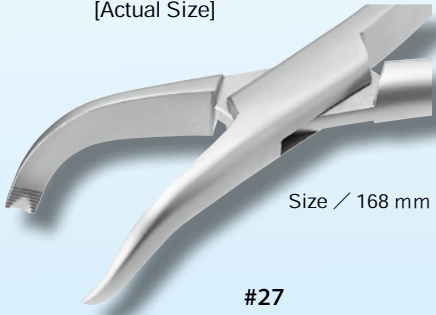

#27
No.17584

Size / 168 mm



[Actual Size]




Po-U
No.17592

Size / 146 mm




[Actual Size]

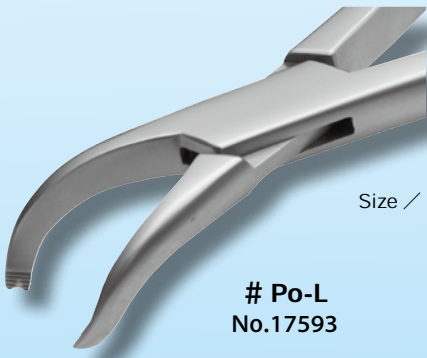

Po-L
No.17593

Size / 146 mm

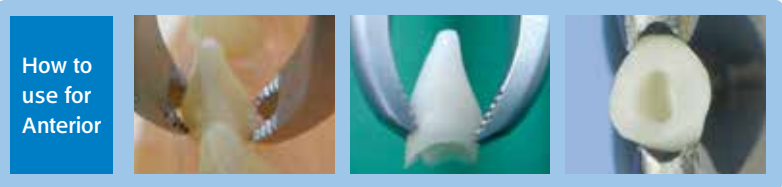

[Actual Size]




Po-U
No.17592

Size / 146 mm



[Actual Size]

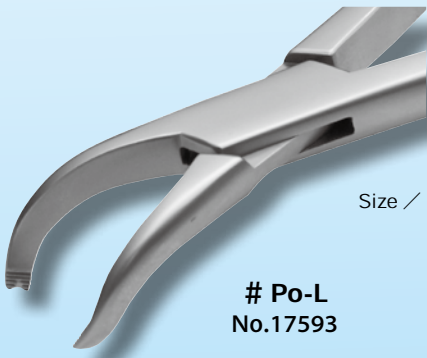

Po-L
No.17593

Size / 146 mm

